# FreeStyle *Libre* 2

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## **App Symbols**

## **Important Safety Information**

Indications for Use

Contraindication

Warnings

Interfering Substances

## **Additional System Information**

## **FreeStyle Libre 2 App Overview**

Home Screen

Sensor Kit

## **App Setup**

## **Applying Your Sensor**

## **Starting Your Sensor**

## **Checking Your Glucose**

## **Understanding Your Glucose Readings**

## **Making Treatment Decisions**

## **Alarms**

## **Setting Alarms**

## **Using Alarms**

## **Adding Notes**

## **Reviewing Your History**

Logbook

Insulin Pens

Other History Options

## **Removing Your Sensor**

## **Replacing Your Sensor**

## **Setting Reminders**

## **Settings and Other Options in the Main Menu**

## **Activities**

## **Maintenance and Disposal**

## **Troubleshooting**

Problems at the Sensor Application Site

Problems Starting Your Sensor or Receiving Sensor Readings

Problems Receiving Glucose Alarms

## **Customer Service**

## **Labelling Symbols and Definitions**

## **Electromagnetic Compatibility**

## **Performance Characteristics**

## User's Manual

# App Symbols


App icon

---

Alarms you have turned on are unavailable

---

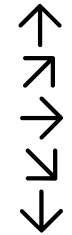Direction your glucose is going. See [Understanding Your Glucose Readings](#) for more information.

---

Caution

---

Add/edit notes

---

Food note

---

Insulin note

---

Food + insulin note

---

Insulin pen error

---

Exercise note

---

Time change

---

Multiple/Custom notes

---

Share report

---

Additional information

---

Main menu

---

Calendar

---

Sensor too cold

---

Sensor too hot

---

## Important Safety Information

### Indications for Use

The FreeStyle Libre 2 Flash Glucose Monitoring System (FreeStyle Libre 2 Reader or FreeStyle Libre 2 app or FreeStyle Libre 2 special edition app used with FreeStyle Libre 2 Sensor) is indicated for measuring interstitial fluid glucose levels in people aged 4 years and older with diabetes mellitus. The System is designed to replace blood glucose testing in the self-management of diabetes, including dosing of insulin. Treatment decisions should be based on real-time Sensor glucose readings and other information on the results screen, including trend arrow and recent sequential readings on the glucose graph. The System also detects trends and tracks patterns and aids in the detection of episodes of hyperglycaemia and hypoglycaemia, facilitating both acute and long-term therapy adjustments.

### Contraindication

**MRI/X-ray/CT:** Remove the Sensor you are wearing if you have a medical appointment that includes strong magnetic or electromagnetic radiation e.g. an X-ray, MRI (Magnetic Resonance Imaging) or CT (Computed Tomography) scan. Apply a new Sensor after the appointment.

## Warnings

- Do not ignore symptoms that may be due to low or high blood glucose. If you have symptoms that do not match the Sensor glucose reading or suspect that your reading may be inaccurate, check the reading by conducting a fingerstick test using a blood glucose meter. If you have symptoms that do not match your glucose readings, consult your health care professional.
- If you are using FreeStyle Libre 2 app ('App'), you must have access to a blood glucose monitoring system as the App does not provide one.
- **Choking hazard:** The System contains small parts that may be dangerous if swallowed.

## Cautions and Important Information

### What to know about Glucose Alarms:

- For you to receive alarms, they must be on and your phone should be within 6 metres (20 ft) of you at all times. The transmission range is 6 metres (20 ft) unobstructed. If you are out of range, you may not receive glucose alarms.
- Do not force close the App.
- Check to make sure that you have the correct settings and permissions enabled on your phone to receive alarms.
  - Enable Bluetooth and allow Nearby Devices permission for the App (Nearby Devices permission is required for Android 12 and above).
  - Allow Notifications for the App. Turn on Lock Screen notifications, Notification sounds, Channel notifications, Sound and Pop-up notifications, and general phone sounds or vibration. Make sure not to turn on any features or modify phone settings that may interrupt the presentation of notifications.
  - Turn off Do Not Disturb mode or manage the settings for your alarms to Override Do Not Disturb. Override Do Not Disturb functionality depends on the phone model and Android operating system version you are using.
  - Allow Alarms & reminders permission for the App. This permission is required for Android 12 and above.
  - Turn off Battery Optimisation for the App to allow the App to run in the background and ensure that you receive alarms, even when your phone has low battery.
  - You may need to add FreeStyle Libre 2 app to the list of apps that will not be put to sleep.

- Remember that alarm settings will follow your phone's sound and vibration settings so these should be at a level you can hear to prevent missed alarms.
- You should disconnect headphones when you are not using them as you may not receive audio with alarms.
- If you are using peripheral devices connected to your phone, such as wireless headphones or a smartwatch, you may receive alarms on only one device or peripheral, not all.
- Keep your phone well charged and turned on.

### **Who should not use the System:**

- The System has not been evaluated for use with other implanted medical devices, such as pacemakers.
- The System has not been evaluated for use in pregnant women, persons on dialysis or people younger than 4 years of age.
- Do not use the System if you are critically ill. The System has not been evaluated for use by critically ill patients and it is not known how different conditions or medications common to this population may affect performance of the System.

### **What to know about wearing the Sensor:**

- Some individuals may be sensitive to the adhesive that keeps the Sensor attached to the skin. If you notice significant skin irritation around or under your Sensor, remove the Sensor and stop using the System. Do not try to treat any skin irritation with lotions/creams/ointments/sprays/barrier patches in order to continue wearing the Sensor. Contact your health care professional before continuing to use the System.
- Intense exercise may cause your Sensor to loosen due to sweat or movement of the Sensor. If your Sensor comes loose, you may get no readings or unreliable readings, which may not match how you feel. Follow the instructions to select an appropriate application site.
- Glucose Going Low and Glucose Going High messages may mean that your actual blood glucose is already  $<3.9$  mmol/L or  $>13.3$  mmol/L. Consider this before making an activity or treatment decision.
- If you are going to participate in activities that could result in harm to yourself or others in the event of a severe hypoglycaemic episode (e.g. driving a motor vehicle without following the Canadian Diabetes Association Recommendations for Private and Commercial Drivers), do not rely on Glucose Messages (High Glucose, Low Glucose, Glucose Going High and Glucose Going Low) alone. Talk to your health care professional

about checking your glucose before and while driving.

- When using the System to replace blood glucose testing for making diabetes treatment decisions, including insulin dosing, you should have a good understanding of how to self-manage your diabetes, as determined by your health care professional.

#### **How to store the Sensor Kit:**

- Store the Sensor Kit between 4°C and 25°C. While you don't need to keep your Sensor Kit in a refrigerator, you can as long as the refrigerator is between 4°C and 25°C.

#### **When not to use the System:**

- Do NOT use if the Sensor Pack or the Sensor Applicator seem to be damaged or already opened.
- Do NOT use if past expiry date.

#### **What to know before you apply the Sensor:**

- Clean the application site and ensure that it is dry prior to Sensor insertion. This helps the Sensor stay attached to your body.

#### **When is Sensor Glucose different from Blood Glucose:**

- Glucose levels in the interstitial fluid can be different from blood glucose levels and may mean that Sensor glucose readings are different from blood glucose. You may notice this difference during times when your blood glucose is changing quickly, for example after eating, taking insulin or exercising.

#### **When to remove the Sensor:**

- If the Sensor is becoming loose or if the Sensor tip is coming out of your skin, you may get no readings or unreliable readings, which may not match how you feel. Check to make sure your Sensor has not come loose. If it has come loose, remove it and apply a new one.
- On rare occasions, you may get inaccurate Sensor glucose readings. If you believe your glucose readings are not correct or do not match how you feel, perform a blood glucose test on your finger to confirm your glucose. If the problem continues, remove your Sensor and apply a new one.

#### **What to know about FreeStyle Libre 2 app:**

- FreeStyle Libre 2 app installed on a phone is intended for use by a single person. It must not be used by more than one person due to the risk of misinterpreting glucose information.

## Interfering Substances

Taking more than 500 mg of ascorbic acid per day may affect the Sensor readings which could cause you to miss a low glucose event. Ascorbic acid can be found in supplements including multivitamins. Some supplements, including cold remedies, may contain high doses of 1000 mg of ascorbic acid. Do not use the Sensor while taking supplements that contain high doses of 1000 mg of ascorbic acid.

## Additional System Information

### App:

- The FreeStyle Libre 2 Reader is a dedicated medical device. Even though you can run FreeStyle Libre 2 app on your phone, the phone is not a dedicated medical device and therefore the App may miss an alarm - for example, due to phone settings, phone or App shutting off, low battery, etc.
- You are responsible for properly securing and managing your phone. If you suspect an adverse cybersecurity event related to FreeStyle Libre 2 app, contact Customer Service.
- Make sure that your phone and Sensor kit are kept in a safe place, under your control. This is important to help prevent anyone from accessing or tampering with the System.
- FreeStyle Libre 2 app is not intended for use on a phone that has been altered or customised to remove, replace or circumvent the manufacturer's approved configuration or use restriction, or that otherwise violates the manufacturer's warranty.

### Sensor:

- The Sensor Pack and Sensor Applicator are packaged as a set and have the same Sensor code. Check that the Sensor codes match before using your Sensor Pack and Sensor Applicator. Sensor Packs and Sensor Applicators with the same Sensor code should be used together or your Sensor glucose readings may be incorrect.
- Your Sensor has been tested to withstand immersion into 1 metre (3 ft) of water for up to 30 minutes. It is also protected against insertion of objects >12 mm diameter. (IP27)
- Do not freeze the Sensor. Do not use if expiry date has passed.
- The Sensor Pack is sterile unless opened or damaged.
- Do not reuse Sensors. The Sensor has been designed so that it cannot be reused. Not suitable for re-sterilisation.

## FreeStyle Libre 2 App Overview

**IMPORTANT:**


- Before you use FreeStyle Libre 2 app with a Sensor, review all the product instructions. The Quick Reference Guide gives you quick access to important aspects of the System. The User's Manual includes all safety information and instructions for use. Refer to your phone's instructions for use for how to use your phone.
- Go to [www.FreeStyleLibre.com](http://www.FreeStyleLibre.com) to view the 'Tips for Kids'.
- Talk to your health care professional about how you should use your Sensor glucose information to help manage your diabetes.
- FreeStyle Libre 2 app is available for download from the Google Play Store. When you're ready to start using FreeStyle Libre 2 app, you'll prepare and apply a Sensor to the back of your upper arm. The App only works with FreeStyle Libre 2 Sensors and cannot be used with other Sensors. When they're in range, the Sensor and App automatically communicate to give you glucose alarms. These alarms are on by default.

**Note:**

- The Sensor comes in a FreeStyle Libre 2 Sensor kit. See [Sensor Kit](#).
- The App is not compatible with all phones. Before upgrading your phone or its operating system, check [www.FreeStyleLibre.com](http://www.FreeStyleLibre.com).

## Home Screen

The Home Screen gives you access to information about your glucose and the App. To return to the Home Screen from another screen, go to the Main Menu and tap **Home**.

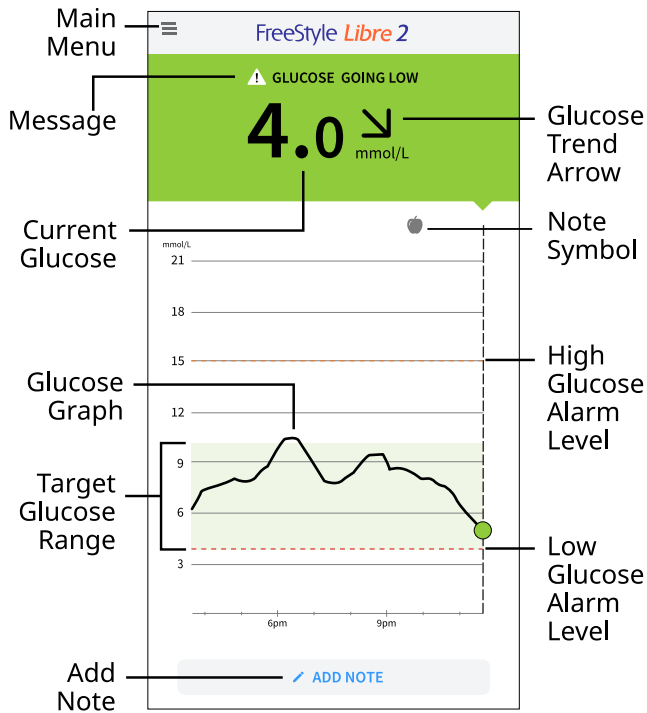**Main Menu** - Tap to access the Home Screen, Alarms, Logbook, other history options and Connected Apps. You can also access Settings, Help and other information.

**Glucose Graph** - Graph of your Sensor glucose readings.

**Current Glucose** - Your most recent glucose value.

**Glucose Trend Arrow** - Direction your glucose is going.

**Message** - You may be able to tap the message for more information.

**Target Glucose Range** - The graph shows your target glucose range. This is not related to glucose alarm levels.

**High Glucose Alarm Level** - Your High Glucose Alarm level displays only if you have kept the alarm on.

**Low Glucose Alarm Level** - Your Low Glucose Alarm level displays only if you have kept the alarm on.

**Add Note** - Tap to add notes to the glucose reading.

**Note Symbol** - Tap to review notes that you've entered.

## Sensor Kit

The Sensor Kit includes:


- Sensor Pack
- Sensor Applicator
- Alcohol wipe
- Product insert

When opening your kit, check that the contents are undamaged and that you have all parts listed. If any parts are missing or damaged, contact Customer Service. The Sensor (only visible after applied) is initially in two parts: one part is in the Sensor Pack and the other part is in the Sensor Applicator. Once prepared and applied to your body, the Sensor measures your glucose using a small, flexible tip that inserts just under the skin. The Sensor can be worn for up to 14 days.

**Sensor Pack.** Used with the Sensor Applicator to prepare the Sensor for use.

**Sensor Applicator.** Applies the Sensor to your body.
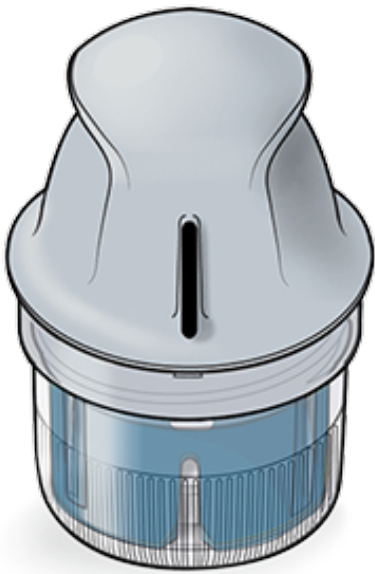
## App Setup

Before using the App for the first time, you must complete the setup.

1. Check that your phone is connected to a network (Wi-Fi or mobile). You can then install FreeStyle Libre 2 from the Google Play Store. Tap the App icon to open the App.  
**Note:** You only need to be connected to a network for setup, using LibreView and sharing with other apps. You do not need to be connected to check your glucose, add notes or review your history in the App.
2. Swipe to view some helpful tips or tap **GET STARTED NOW** at any point.

3. Confirm your country/region and tap **NEXT**.
4. You need a LibreView account to use the App. Follow onscreen instructions to review legal information and create a new account or login to your existing account.
5. Confirm your glucose unit of measure and tap **NEXT**.
6. Select how you count carbohydrates (in grams or portions) and tap **NEXT**. The carbohydrate unit will be used in any food notes that you enter in the App.
7. Select whether you want sound and vibration OR vibration only when you scan your Sensor. Tap **NEXT**.

**Note:** This setting does not affect alarms.


8. The App now displays some useful information. Tap **NEXT** to review each screen.
9. Accept required permissions.
10. Apply a new Sensor and then tap **NEXT**. Go to [Starting Your Sensor](#).

**Note:** If you need help applying your Sensor, tap **HOW TO APPLY A SENSOR** or go to [Applying Your Sensor](#).

## Applying Your Sensor

### CAUTION:

- The Sensor Pack and Sensor Applicator are packaged together as a set and have the same Sensor code. Check that the Sensor codes match before using your Sensor Pack and Sensor Applicator. Sensor Packs and Sensor Applicators with the same Sensor code should be used together or your Sensor glucose readings may be incorrect.

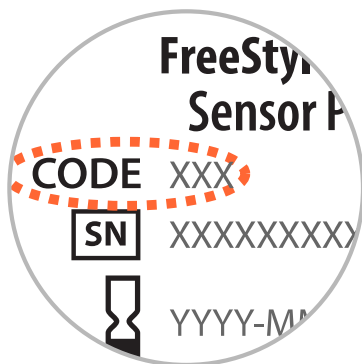- Intense exercise may cause your Sensor to loosen due to sweat or movement of the Sensor. If your Sensor comes loose, you may get no readings or unreliable readings, which may not match how you feel. Follow the instructions to select an appropriate

application site.


---

1. Apply Sensors only on the back of your upper arm. Avoid areas with scars, moles, stretch marks or lumps. Select an area of skin that generally stays flat during your normal daily activities (no bending or folding). Choose a site that is at least 2.5 cm (1 inch) away from an insulin injection site. To prevent discomfort or skin irritation, you should select a different site other than the one most recently used.


2. Wash application site for at least 10 seconds using a plain soap and dry with a towel. Make sure site is dry and then clean with an alcohol wipe. This will help remove any oily residue that may prevent the Sensor from sticking properly. Allow site to air dry before proceeding.

**Note:** The area **MUST** be clean and dry, or the Sensor may not stick to the site.

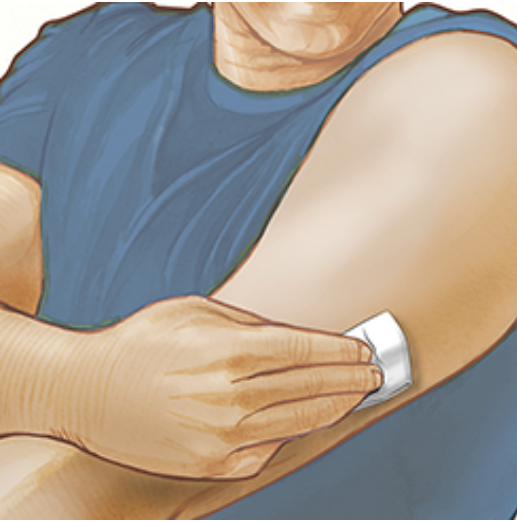
3. Open the Sensor Pack by peeling the lid off completely. Unscrew the cap from the Sensor Applicator and set the cap aside.

**CAUTION:** Do NOT use if the Sensor Pack or Sensor Applicator seem to be damaged or already opened. Do NOT use if past expiry date.

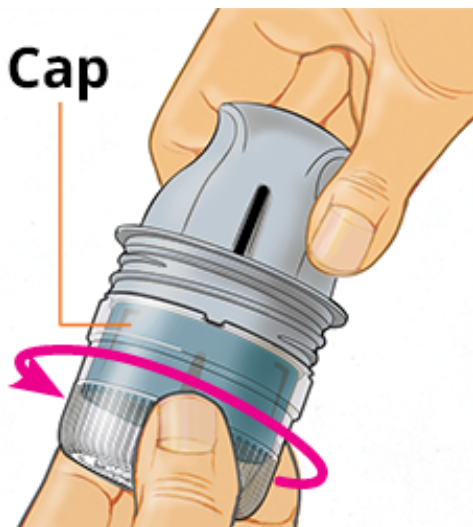4. Line up the dark mark on the Sensor Applicator with the dark mark on the Sensor Pack. On a hard surface, press down firmly on the Sensor Applicator until it comes to a stop.

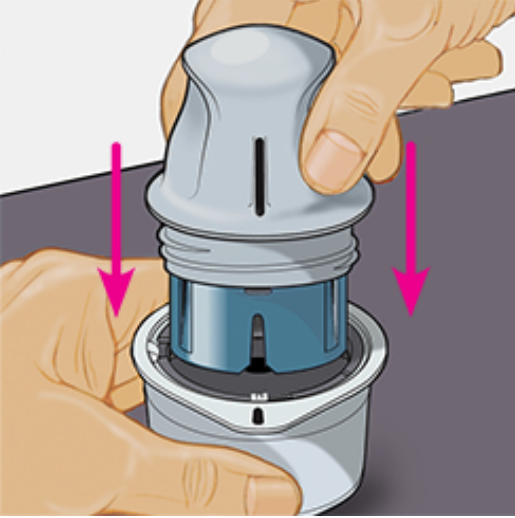
5. Lift the Sensor Applicator out of the Sensor Pack.

6. The Sensor Applicator is prepared and ready to apply the Sensor.


**CAUTION:** The Sensor Applicator now contains a needle. Do NOT touch inside the Sensor Applicator or put it back into the Sensor Pack.

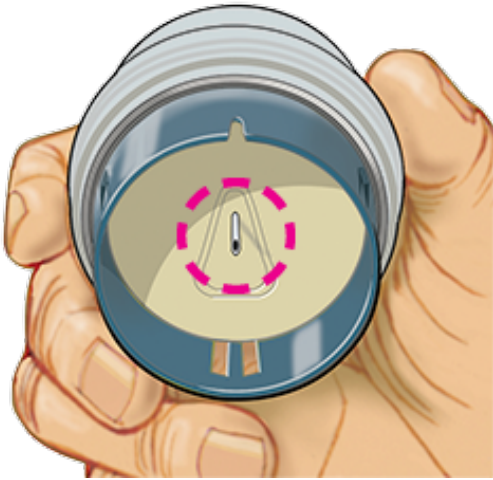
7. Place the Sensor Applicator over the prepared site and push down firmly to apply the Sensor to your body.

**CAUTION:** Do NOT push down on Sensor Applicator until placed over prepared site to prevent unintended results or injury.

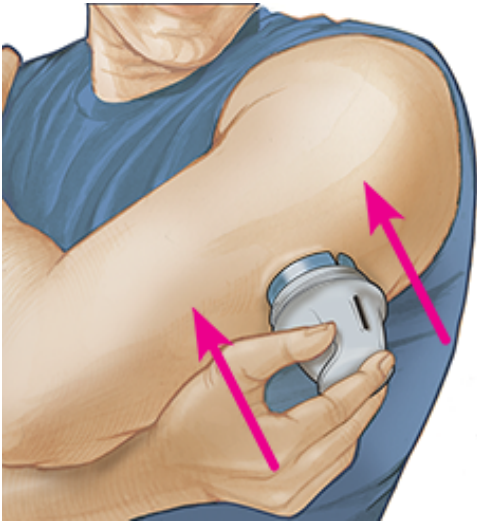8. Gently pull the Sensor Applicator away from your body. The Sensor should now be attached to your skin.

**Note:** Applying the Sensor may cause bruising or bleeding. If there is bleeding that does not stop, remove the Sensor, and apply a new one at a different site.

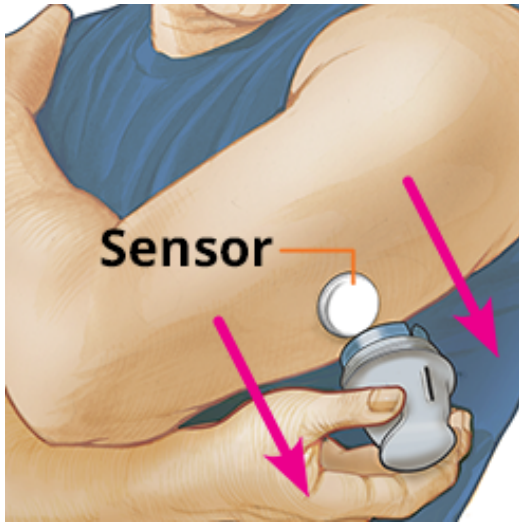9. Make sure Sensor is secure after application. Put the cap back on the Sensor Applicator. Discard the used Sensor Applicator and Sensor Pack according to local regulations.

**Note:** Tap **Help** in the Main Menu to access an in-app tutorial on applying a Sensor.

## Starting Your Sensor

### IMPORTANT:

- Keep in mind that the ease of scanning a Sensor may vary between devices. Once you determine the location of the NFC (Near Field Communication) antenna on your phone, you should be able to reliably scan your Sensor by holding that area near your Sensor. You may need to adjust your scan distance based on what clothing you are wearing. In addition to proximity and orientation, other factors can affect NFC performance. For example, a bulky or metallic case can interfere with the NFC signal.
- The App requires that your phone be set to network time. This is the default time setting for most phones.

- When using the App, you should keep your phone well charged and make sure you have access to a blood glucose meter.
  - Make sure your phone's sound is set so you can hear scan sounds or alarms if you turned those on.
- 

1. Hold the back of your phone near the Sensor (this can be done over clothing) and position the NFC antenna over the Sensor. Do not move the phone until you hear the first tone and/or feel a vibration. This indicates your phone and Sensor have established an NFC connection.

**Note:** Before you start your Sensor, choose which device you want to use. If you start the Sensor with the FreeStyle Libre 2 app, you will be unable to use the Reader or FreeStyle Libre 2 special edition app to check your glucose or receive alarms.

2. Continue holding your phone near the Sensor until you hear a second tone and/or feel a vibration. This completes the scan.

**Note:**

- If you need help, tap **HOW TO SCAN A SENSOR** to view an in-app tutorial. You can also access this later by going to the Main Menu and then tapping **Help**.
- If your Sensor is not successfully scanned, you may receive a scan error.

See [Troubleshooting](#) for additional error messages.

3. The Sensor can be used to check your glucose after 60 minutes. While the Sensor is starting up, you can navigate away from the App. You will see a notification when the Sensor is ready.

**CAUTION:** Low and High Glucose Alarms are not available during the 60-minute start-up period.

**Note:** If you get a new phone after you've started a Sensor, you can continue to use the Sensor and receive alarms. Your settings and data are tied to your LibreView account and not to your phone so you just need to download the App again and log into your LibreView account.


## Checking Your Glucose

1. You'll see your glucose reading whenever you open the App, automatically updated every minute. If your glucose reading does not automatically display, update the App and continue to scan your current Sensor to get glucose readings. When you start your next Sensor with the updated App, you should see your readings display and automatically

update every minute.

You can still scan your Sensor whenever you want – for example, to fill in up to 8 hours of missing data or to get a glucose reading during signal loss. To get your glucose reading with a scan, just hold the back of your phone near the Sensor. If scan sounds are turned on, you will get two separate tones along with the vibrations when your Sensor has been scanned.

- Your glucose reading includes your Current Glucose, a Glucose Trend Arrow indicating which way your glucose is going and a graph of your current and stored glucose readings.

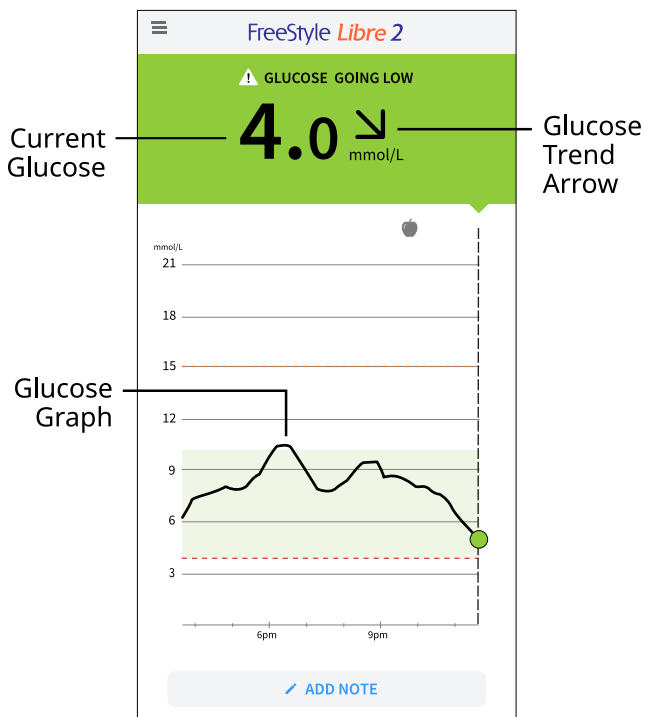**Current Glucose** - Your most recent glucose value.

**Glucose Trend Arrow** - Direction your glucose is going.

**Glucose Graph** - Graph of your Sensor glucose readings.

#### Note:

- A Sensor can store up to 8 hours of glucose data. Scan it when you see gaps in your graph to capture all of your available glucose data.
- The graph will scale to 22.2 mmol/L to accommodate glucose readings above 21 mmol/L.
- The ⌚ symbol may appear, indicating your phone's time was changed. Gaps in the graph may result or glucose readings may be hidden.

- All available glucose data is used to make your graph so you can expect to see some differences between the graph line and previous current glucose readings.
- Your current glucose value determines the background colour on your glucose reading display:

**Orange** - High glucose (above 13.3 mmol/L)

**Yellow** - Between the Target Glucose Range and high or low glucose level

**Green** - Within the Target Glucose Range

**Red** - Low glucose (below 3.9 mmol/L)

## Understanding Your Glucose Readings

### Glucose Trend Arrow

The Glucose Trend Arrow gives you an indication of the direction your glucose is going.

Glucose is rising quickly (more than 0.1 mmol/L per minute)

Glucose is rising (between 0.06 and 0.1 mmol/L per minute)

Glucose is changing slowly (less than 0.06 mmol/L per minute)

Glucose is falling (between 0.06 and 0.1 mmol/L per minute)


Glucose is falling quickly (more than 0.1 mmol/L per minute)


### Messages

Below are messages you may see with your glucose readings.

**LO | HI:** If **LO** appears, your reading is lower than 2.2 mmol/L. If **HI** appears, your reading is higher than 22.2 mmol/L. You can tap the  symbol for more information. Check your blood glucose on your finger with a test strip. If you get a second **LO** or **HI** result, contact your health care professional IMMEDIATELY.

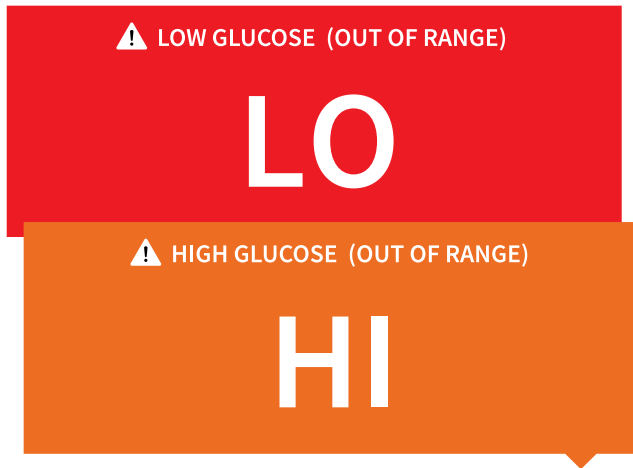**Low Glucose | High Glucose:** If your glucose is higher than 13.3 mmol/L or lower than 3.9 mmol/L, you will see a message on the screen. You can tap the  symbol for more information and set a reminder to check your glucose.

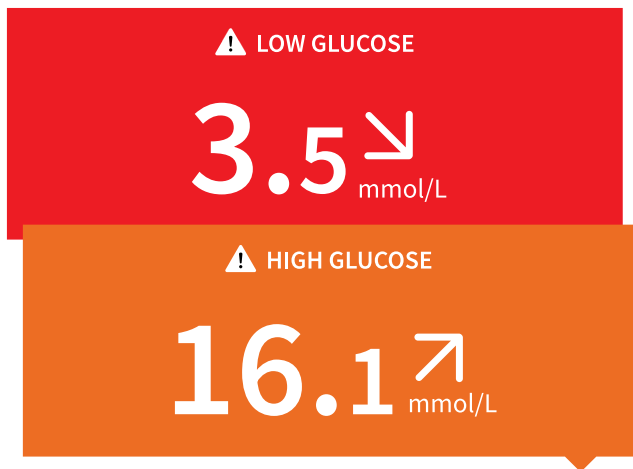**Glucose Going Low | Glucose Going High:** If your glucose is projected to be higher than 13.3 mmol/L or lower than 3.9 mmol/L within 15 minutes, you will see a message on the screen. The background colour corresponds to your current glucose value. You can tap the  symbol for more information and set a reminder to check your glucose.

A green rectangular box with a white warning triangle icon on the left. The text "GLUCOSE GOING LOW" is in white. Below it, "4.0" is written in large black font with a downward-pointing arrow to its right. Underneath "4.0" is "mmol/L" in a smaller black font.

⚠ GLUCOSE GOING LOW

4.0  $\searrow$   
mmol/L

An orange rectangular box with a white warning triangle icon on the left. The text "GLUCOSE GOING HIGH" is in white. Below it, "13.2" is written in large black font with an upward-pointing arrow to its right. Underneath "13.2" is "mmol/L" in a smaller black font.

⚠ GLUCOSE GOING HIGH

13.2  $\nearrow$   
mmol/L

### Note:

- If you are not sure about a message or reading, contact your health care professional for information.
- Messages you receive with glucose readings (**LO, HI, Low Glucose, High Glucose, Glucose Going Low, Glucose Going High**) are not related to glucose alarm settings.

## Making Treatment Decisions


Work with your health care professional to put together a plan for managing your diabetes that includes when to use the System information for making treatment decisions.

### WARNING:

The System can replace blood glucose testing except in the below situations. These are the times when you need to do a blood glucose test before deciding what to do or what treatment decision to make as Sensor readings may not accurately reflect blood glucose levels:

**Do a blood glucose test** if you think your glucose readings are not correct or do not match how you feel. Do not ignore symptoms that may be due to low or high glucose.

**Do a blood glucose test** if the Sensor glucose reading does not include a Current Glucose number.
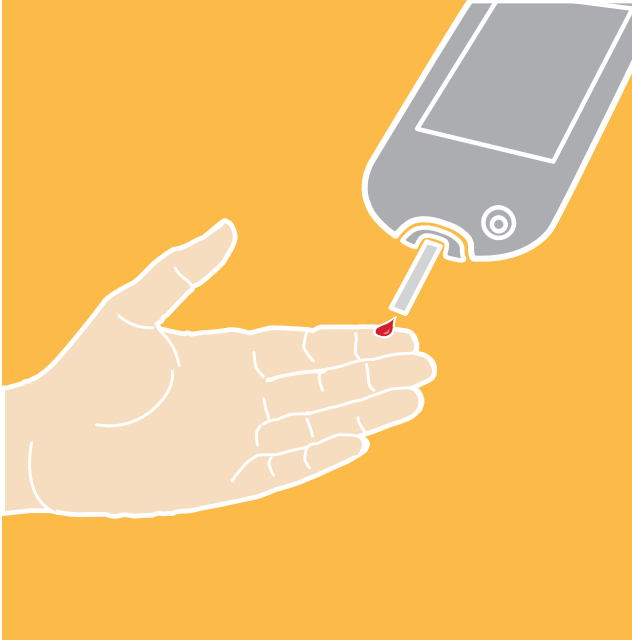
## Making Treatment Decisions – Getting Started

Before you start using the System for treatment decisions, make sure that you have a good understanding of how the System works for your body. **Continue to use your blood glucose meter for treatment decisions until you are comfortable with the information you receive from your System.** This includes understanding that: Sensor performance can vary in between Sensors, within a Sensor wear period and in different situations.

Getting familiar with the System could take days, weeks or even months. The more you check readings from the System with a blood glucose meter, the better you will understand how the System works for you.

Work with your health care professional to put together a plan for managing your diabetes that includes when to use the System information for making treatment decisions.

### Helpful Tips

- Confirm your Sensor glucose readings with a blood glucose meter until you understand:
  - Sensor accuracy may vary between Sensors. Some individual Sensors may read higher or lower than true blood glucose levels. This may be temporary. For example, when glucose is going up or down quickly. Alternatively, it could apply to all Sensor glucose readings over the 14-day wear period due to variability in factory calibration and changes in the Sensor during storage.
  - Sensor accuracy may vary during a Sensor wear session.
  - Sensor accuracy may vary in different situations (meals, exercise, first day of use, etc.). Sensors may tend to give lower and more variable glucose readings on day 1 of Sensor

wear since it takes some time for your body to adjust to a newly inserted Sensor.

- Check your glucose often to see how carbs, medication, exercise, illness or stress levels impact your Sensor glucose readings. The information you get can help you figure out why your glucose sometimes goes too high or too low, and how to prevent it from doing so in the future.
- Talk to your health care professional about how your insulin works. The more you understand about your insulin, including how long it takes to start working and how long it lasts in your body, the more likely you will be to make better treatment decisions.
- Making a treatment decision doesn't just mean taking insulin. Treatment decisions can also include things like taking fast-acting carbs, eating or even doing nothing and checking again later.
- Your health care professional can also help you to understand when doing nothing and checking your glucose again later is the right treatment decision. For example, if your glucose is high and going up, your first instinct may be to take more insulin to lower your glucose. However, depending on when you last took insulin or your recent activity, the right treatment decision may be to do nothing and check again later. Avoid 'insulin stacking'.

## When not to use Sensor Glucose readings for treatment decisions

### ***No Current Glucose Number***

When there is no Current Glucose number, such as when you receive an error message or a LO or HI result, you don't have enough information to make a treatment decision. Do a blood glucose test and treat based on that result.

### ***Think Your Readings are Incorrect?***

Don't trust Sensor glucose readings that you think may be incorrect or that don't match what you would expect based on your recent activity. For example, if you ate dinner but forgot to take insulin before eating, you would expect your glucose to be high. If your glucose reading is low, then it doesn't match your recent activity, so don't use it to make treatment decisions. Don't make treatment decisions if you think your Sensor glucose readings are incorrect. Do a blood glucose test and treat based on that result.

### ***Symptoms Don't Match Readings***

There may be times when your symptoms don't match your Sensor glucose readings. For example, you are feeling shaky, sweaty and dizzy – symptoms you generally get when you have low glucose, but your glucose reading is within your target range. When symptoms don't match readings, do a blood glucose test and treat based on that result. Don't ignore

symptoms that may be due to low or high blood glucose.

If you're the carer, pay attention to times when the symptoms of the one you're caring for don't match their Sensor glucose readings. When symptoms don't match readings, do a blood glucose test and treat based on that result.

If you continue to feel that Sensor glucose readings are not correct, compare the blood glucose result to the Sensor glucose reading. You'll know the Sensor is performing correctly if these results are within 20% of each other (for example, when the Sensor reads 10.0 mmol/L, the blood glucose result should be between 8.0 and 12.0 mmol/L). If the blood glucose result is not within 20% of the Sensor glucose reading, use a blood glucose meter instead of the Sensor and try the check again later. If the problem continues, keep using the blood glucose meter instead of the Sensor and contact your health care professional for guidance on when you can use the Sensor.


### ***When to do Nothing and Check Your Glucose Again Later***

Your health care professional can help you understand when doing nothing and checking your glucose again later is the right treatment decision. For example, if your glucose is high and going up, your first instinct may be to take more insulin to lower your glucose. However, depending on when you last took insulin or your recent activity, the right treatment decision may be to do nothing and check your glucose again later.

Don't take a correction dose within 2 hours of your meal dose. This may result in 'insulin stacking' and low glucose.

## **Understanding the Difference between ISF Glucose and Blood Glucose**

Please contact your health care professional if you need help understanding the information below. The Sensor measures interstitial fluid (ISF) instead of blood glucose. ISF glucose levels tend to lag behind blood glucose levels because it can take a few minutes for glucose to move from the blood into the ISF. This is particularly true when blood glucose levels are changing quickly. As a result, when glucose levels are falling quickly, glucose readings from the Sensor may be higher than blood glucose levels. Conversely, when glucose levels are rising quickly, glucose readings from the Sensor may be lower than blood glucose levels.
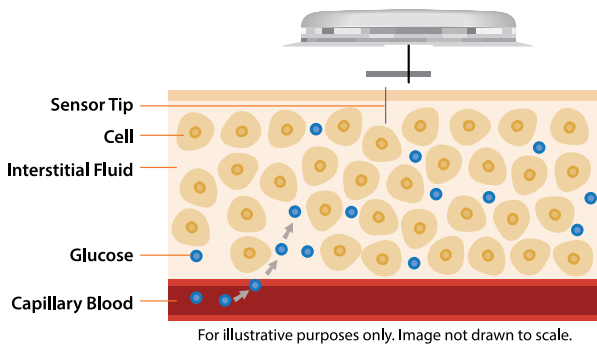
## Using Your Glucose Reading to Make a Treatment Decision

After you check your glucose, use all of the information on the screen when deciding what to do or what treatment decision to make.

**Glucose Trend Arrow**  
Direction your glucose is going

| Arrow | What it means           |
|-------|-------------------------|
| ↑     | Glucose rising quickly  |
| ↗     | Glucose rising          |
| →     | Glucose changing slowly |
| ↘     | Glucose falling         |
| ↓     | Glucose falling quickly |

**Message**  
Tap for more information


**Current Glucose**  
Your most recent glucose value

**Glucose Graph**  
Graph of your current and stored glucose readings

### Example Scenarios

Here are some example scenarios to help you understand how to use the information on your screen. Always use all of the information on the screen before deciding what to do or treatment decision to make. If you are not sure about what to do, consult your health care professional.

## What you see – *Before breakfast*
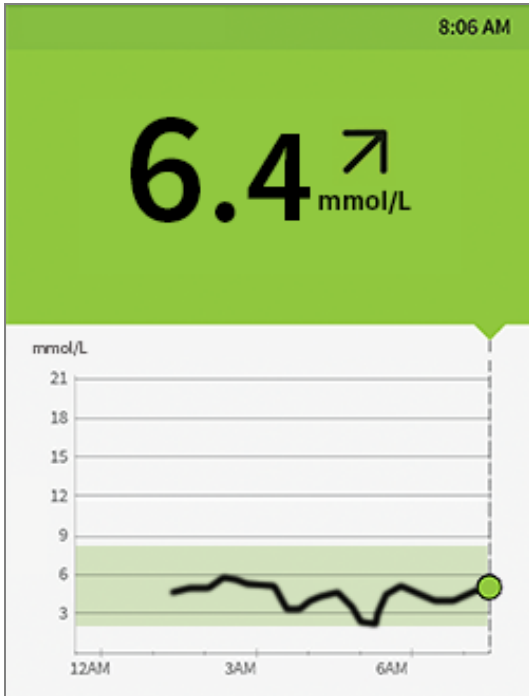
## What it means

Before breakfast, your current glucose is 6.4 mmol/L. The graph shows that your glucose is going up and so does the trend arrow ↗.

Consider what might be causing your glucose to go up and what you might do to prevent a high glucose. For example:

- How much insulin should you take before your meal?
- Since you see ↗, should you consider taking a little more insulin?

## What you see – *Before lunch*

### What it means

When you checked your glucose before lunch, it was 5.0 mmol/L and rising. Before eating lunch, you took enough insulin to cover the meal and a little more since your trend arrow was  $\nearrow$ .

### What you see - After lunch

### What it means

90 minutes later, your current glucose is 12.5 mmol/L. The graph shows that your glucose is still going up, and so does the trend arrow  $\nearrow$ .

**Don't take a correction dose within 2 hours of your meal dose.** This may result in 'insulin stacking' and low glucose.

Consider what might be causing your glucose to go up and what you might do to prevent a high glucose. For example:

- Has the insulin you took for your meal reached its full effect?
- Check your glucose again later.

### What you see - *In the afternoon*


### What it means

Between meals, your current glucose is 4.0 mmol/L. The Glucose Going Low message tells you that your glucose is projected to be low within 15 minutes.

Think about what might be causing your glucose to go low. Consider eating a snack to stay within target. **Avoid taking insulin as this can cause low glucose.**

### What you see - *After exercising*
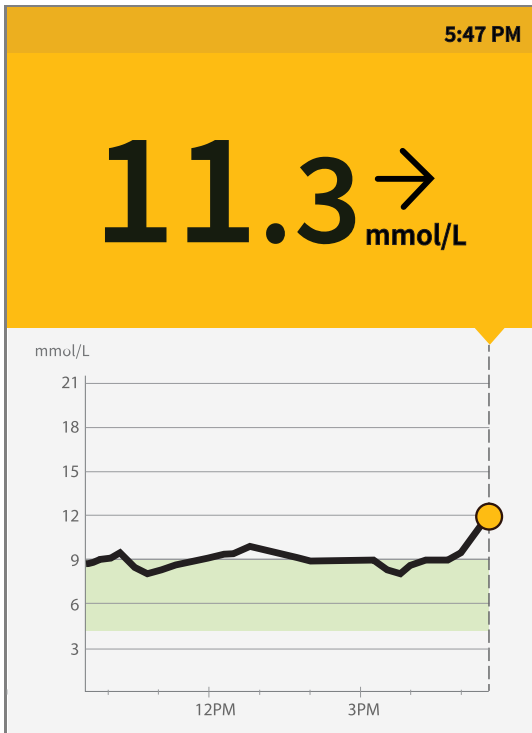
### What it means

After exercising, you are feeling shaky, sweaty and dizzy – symptoms you generally get when you have low glucose. But, your current glucose is 11.3 mmol/L.

Any time that you get a reading that doesn't match how you feel, do a blood glucose test.

### What you see – *Before dinner*

### What it means

Before dinner, your current glucose is 7.4 mmol/L. The graph shows that your glucose is going down and so does the trend arrow ↘.

Consider what might be causing your glucose to go down and what you might do to prevent a low glucose. For example:

- How much insulin should you take to cover your meal?
- Since you see ↘, should you think about taking a little less insulin?

## Alarms

When in range of your phone, your Sensor automatically communicates with the App to give you Low and High Glucose Alarms. These alarms are on by default.

This section explains how to set and use alarms as well as how to turn them off. Please read all the information in this section before setting and using alarms.

### CAUTION:

- For you to receive alarms, they must be on and your phone should be within 6 metres (20 ft) of you at all times. The transmission range is 6 metres (20 ft) unobstructed. If you are out of range, you may not receive glucose alarms.
- Do not force close the App.
- Check to make sure that you have the correct settings and permissions enabled on your phone to receive alarms.
  - Enable Bluetooth and allow Nearby Devices permission for the App (Nearby Devices permission is required for Android 12 and above).
  - Allow Notifications for the App. Turn on Lock Screen notifications, Notification sounds, Channel notifications, Sound and Pop-up notifications, and general phone sounds or vibration. Be sure not to turn on any features or modify phone settings that may interrupt the presentation of notifications.
  - Turn off Do Not Disturb mode or manage the settings for your alarms to override Do Not Disturb. Override Do Not Disturb functionality depends on the phone model and Android operating system version you are using.
  - Allow Alarms & reminders permission for the App. This permission is required for Android 12 and above.
  - Turn off Battery Optimisation for the App to allow the App to run in the background and ensure that you receive alarms, even when your phone has low battery.

- You may need to add FreeStyle Libre 2 app to the list of apps that will not be put to sleep.
- Remember that alarm settings will follow your phone's sounds and vibration settings so these should be at a level you can hear to prevent missed alarms.
- You should disconnect headphones when you are not using them as you may not receive audio with alarms.
- If you are using peripheral devices connected to your phone, such as wireless headphones or a smartwatch, you may receive alarms on only one device or peripheral, not all.
- Keep your phone well charged and turned on.

### **IMPORTANT: What to know about glucose alarms**

- The Sensor itself will not issue alarms.
- The Low and High Glucose Alarms should not be used exclusively to detect low or high glucose conditions. The glucose alarms should always be used along with your current glucose, glucose trend arrow and glucose graph.
- Low and High Glucose Alarm levels are different from your Target Glucose Range values. Low and High Glucose Alarms tell you when your glucose has passed the level that you set in the alarm. Your Target Glucose Range is displayed on glucose graphs in the App and used to calculate your Time In Range.
- Make sure that the **Signal Loss Alarm** is turned on so you will be notified if you are not receiving your Low or High Glucose alarms.
- You will not receive glucose alarms if you completely close the App.

### **IMPORTANT: How to prevent missed alarms**

- Make sure that your phone is near you. The Sensor itself will not issue alarms.
- If the Sensor is not communicating with the App, you will not receive glucose alarms, and you may miss detecting low glucose or high glucose episodes. You will see the  or  symbol on the screen when the Sensor is not communicating with the App. The  or  symbol will also display when the System is unable to calculate a glucose reading due to an error. A message on the screen will give you directions to resolve the error. Make sure that the **Signal Loss Alarm** is on so you will be notified if your Sensor has not communicated with the App for 20 minutes.
- If you see the  or  symbol, this means that you are not getting glucose alarms because of one or more of the following:


- Bluetooth is OFF
- App notifications are OFF
- Nearby Devices permission for the App is OFF. This permission is required for Android 12 and above.
- Alarms & reminders permission for the App is OFF. This permission is required for Android 12 and above.
- Sensor is not communicating with the App or glucose reading error
- Lock Screen notifications, or Notification sounds are OFF
- Channel notifications or Sound and Pop-up notifications are OFF
- Battery optimisation is ON

## Setting Alarms

To set or turn off alarms, go to the Main Menu and tap **Alarms**. Select the alarm you want to change.


### Low Glucose Alarm

1. The **Low Glucose Alarm** is on by default. The alarm level is initially set to 3.9 mmol/L. Tap to change this value between 3.3 mmol/L and 5.6 mmol/L. If the alarm is on, you will be notified when your glucose falls below the level you set. Touch the slider to turn the alarm off. Tap **SAVE**.
2. Choose the sound for this alarm. Volume and vibration will match your phone settings.
3. Select whether to manage the settings for this alarm to **Override Do Not Disturb**. Override Do Not Disturb functionality depends on the phone model and Android operating system version you are using.
4. Tap the back button to return to the main alarm settings screen.

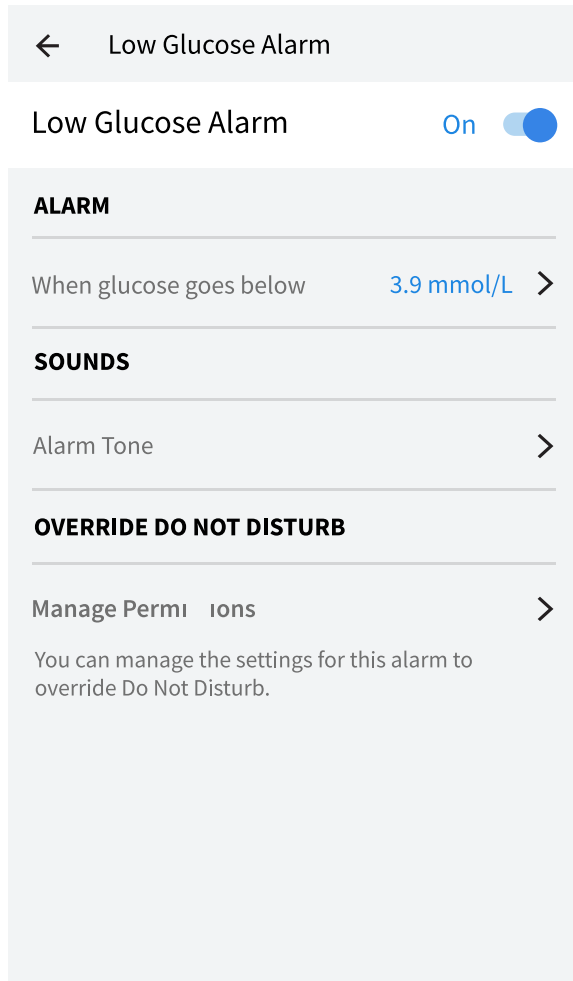
## High Glucose Alarm

1. The **High Glucose Alarm** is on by default. The alarm level is initially set to 13.3 mmol/L. Tap to change this value between 6.7 mmol/L and 22.2 mmol/L. If the alarm is on, you will be notified when your glucose rises above the level you set. Touch the slider to turn the alarm off. Tap **SAVE**.
2. Choose the sound for this alarm. Volume and vibration will match your phone settings.
3. Select whether to manage the settings for this alarm to **Override Do Not Disturb**. Override Do Not Disturb functionality depends on the phone model and Android operating system version you are using.
4. Tap the back button to return to the main alarm settings screen.
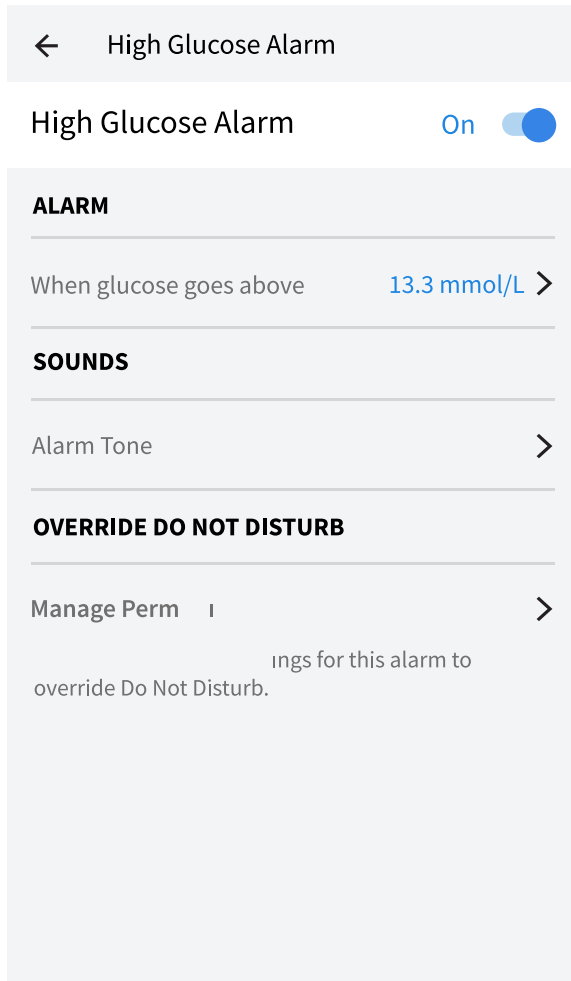
## Signal Loss Alarm

1. If the alarm is on, you will be notified when your Sensor has not communicated with the App for 20 minutes and you are not receiving Low or High Glucose Alarms. Tap the slider to turn the alarm off.
2. Choose the sound for this alarm. Volume and vibration will match your phone settings.
3. Select whether to manage the settings for this alarm to **Override Do Not Disturb**. Override Do Not Disturb functionality depends on the phone model and Android operating system version you are using.
4. Tap the back button to return to the main alarm settings screen.
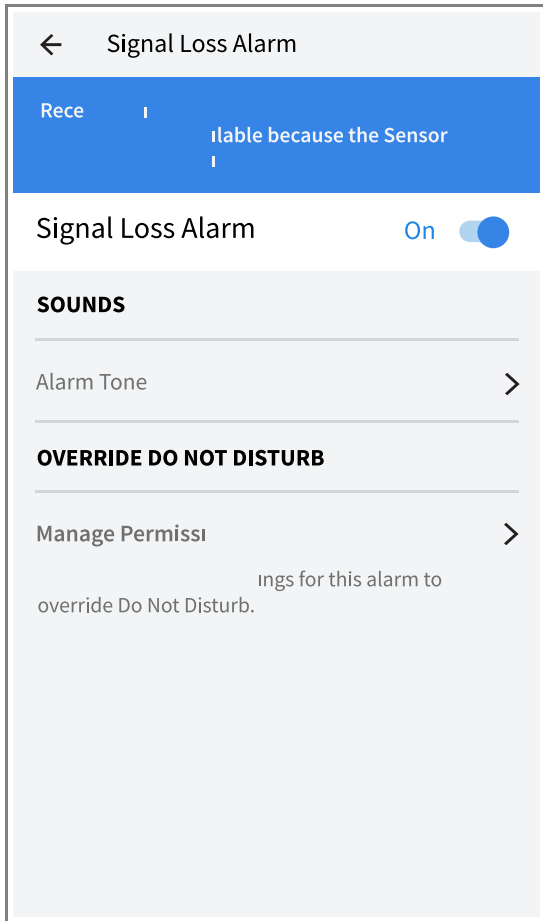
## Using Alarms

**Low Glucose Alarm** notifies you if your glucose drops below the level you set. Open the App or tap the **Dismiss** button to dismiss the alarm. You will only receive one alarm per low glucose episode.

Low Glucose Alarm ⚠  
3.7 mmol/L ↘

**High Glucose Alarm** notifies you if your glucose rises above the level you set. Open the App or tap the **Dismiss** button to dismiss the alarm. You will only receive one alarm per high glucose episode.

High Glucose Alarm ⚠  
13.4 mmol/L ↗

**Signal Loss Alarm** notifies you if your Sensor has not communicated with the App for

20 minutes and you are not receiving Low or High Glucose Alarms. Signal loss could be caused by the Sensor being too far away from your smartphone (over 6 metres (20 ft)) or another issue such as an error or problem with your Sensor. Open the App or tap the **Dismiss** button to dismiss the alarm.

Signal Loss Alarm   
Alarms are unavailable. Scan Sensor.

### Note:

- If you ignore an alarm, you will receive it again in 5 minutes if the condition still exists.
- Only your most recent alarms will display on your screen.

## Adding Notes

Notes can be saved with your glucose readings to help you track food, insulin and exercise. You can also add your own comment.

1. Tap  on your glucose reading screen.
2. Select the tick box next to the notes you would like to add. After you tick the box, you can add more specific information to your note.
  - Food notes: Enter meal type and grams or portion information
  - Insulin notes: Enter the number of units taken
  - Exercise notes: Enter intensity and duration
3. Tap **DONE** to save your note.

Notes you add are shown on your glucose graph and in your Logbook as symbols. You can review a note by tapping its symbol on your glucose graph or by going to the Logbook. See [Reviewing Your History](#) for more information about the Logbook. To edit a note from the glucose graph, tap the symbol and then tap on the information you would like to change. Tap **DONE** when you are finished.

Food

Insulin

Exercise

---

Food + insulin

---

Multiple/Custom notes – indicates different types of notes entered together or notes entered within a short period of time. A numbered badge next to the symbol indicates the number of notes.

---

## Reviewing Your History

Reviewing and understanding your glucose history can be an important tool for improving your glucose control. The App stores about 90 days of information and has several ways to review your past glucose readings and notes. From the Main Menu, tap **Logbook** to view the Logbook or tap on one of the other history options under **Reports**.

**IMPORTANT:** Work with your health care professional to understand your glucose history.

---

### Logbook

The **Logbook** contains entries for each time that you scanned your Sensor as well as notes that you added. If you would like to view a different day, tap the  symbol or use the arrows. To add a note to a Logbook entry, tap on the entry and then tap . Select your note information and tap **DONE**.

To add a note that is independent of a Logbook entry, tap  on the main Logbook screen. Tap  if you want to add a note on a different date.

### Insulin Pens

The **Insulin Pens** option in the Main Menu allows you to connect compatible insulin pens to the App. Once you have connected your insulin pen, you can transfer insulin doses from your insulin pen to the App. Insulin doses can be reviewed in the Logbook. For more information, refer to the Insulin Pen User Guide in the Help section. Review all instructions provided by your insulin pen manufacturer prior to insulin pen use.

### Other History Options

**Daily Patterns:** A graph showing the pattern and variability of your Sensor glucose readings

over a typical day. The thick black line shows the median (midpoint) of your glucose readings. The light blue shading represents the 10th – 90th percentile range of your glucose readings. Dark blue shading represents the 25th – 75th percentile range.

**Note:** This report needs at least 5 days of glucose data.

**Time In Range:** A graph showing the percentage of time your Sensor glucose readings were above, below or within your Target Glucose Range.

**Low Glucose Events:** Information about the number of low glucose events measured by your Sensor. A low glucose event is recorded when your Sensor glucose reading is lower than 3.9 mmol/L for 15 minutes or longer. The total number of events is displayed below the graph. The bar graph displays the low glucose events in different periods of the day.

**Average Glucose:** Information about the average of your Sensor glucose readings. The overall average for the selected time period is displayed below the graph. The average is also shown for different periods of the day. Readings above or below your Target Glucose Range are yellow, orange or red. Readings in range are green.

**Daily Graph:** A graph of your Sensor glucose readings by day. The graph shows your Target Glucose Range and symbols for notes you have entered.

- The graph will scale to 22.2 mmol/L to accommodate glucose readings above 21 mmol/L.
- You might see gaps in the graph during times when you have not scanned at least once in 8 hours or if you lost Bluetooth connection.
- The ⌚ symbol may appear indicating a time change. Gaps in the graph may result or glucose readings may be hidden.

**Estimated A1c:** Your estimated A1c level (also called HbA1c) is based on available Sensor glucose data from the last 90 days. The more data available, the better your estimation will be. However, the estimated level may not match your A1c measured in a laboratory\*. A1c can be used as an indicator of how well your glucose levels have been controlled and may be used to monitor your diabetes treatment regimen.

\* The formula is based on the published reference, which compared average Sensor glucose and laboratory-measured A1c:

$$A1c\% = (\text{Avg } SG_{\text{mmol/L}} + 2.59) / 1.59$$

Reference: Nathan DM, Kuenen J, Borg R, Zheng H, Schoenfeld D, Heine RJ for the A1c-Derived Average Glucose (ADAG). Study Group: Translating the hemoglobin A1c assay into estimated average glucose values. Diabetes Care 2008, 31:1473-8.

**Sensor Usage:** Information about how often you scanned your Sensor or viewed your Sensor glucose readings in the App and how much information has been captured from your Sensor.

## Note:

- Tap the  symbol on any report to share a screenshot of the report.
- Tap the  symbol to view a description of the report.
- From any report screen, swipe left or right to view the next or previous report.
- On all reports except the **Daily Graph** and **Estimated A1c**, you can select to show information about your last 7, 14, 30 or 90 days.

## Removing Your Sensor

1. Pull up the edge of the adhesive that keeps your Sensor attached to your skin. Slowly peel away from your skin in one motion.

**Note:** Any remaining adhesive residue on the skin can be removed with warm soapy water or isopropyl alcohol.

2. Discard the used Sensor according to local regulations. When you are ready to apply a new Sensor, follow the instructions in [Applying Your Sensor](#) and [Starting Your Sensor](#).

## Replacing Your Sensor

Your Sensor automatically stops working after 14 days of wear and must be replaced. You should also replace your Sensor if you notice any irritation or discomfort at the application site or if the App reports a problem with the Sensor currently in use. Taking action early can keep small problems from turning into larger ones.

**CAUTION:** If the Sensor is becoming loose or if the Sensor tip is coming out of your skin, you may get no readings or unreliable readings, which may not match how you feel. Check to make sure your Sensor has not come loose. If it has come loose, remove it and apply a new one.

# Setting Reminders

You can create single or repeating reminders to help you remember things like checking your glucose or taking insulin. There is one default reminder to help you remember to check your glucose. This reminder can be changed or disabled but cannot be deleted.

## Note:

- If you want to receive a sound/vibration with your reminder, ensure that sound/vibration on your smartphone is turned on, sound is set at a level you can hear and your smartphone's Do Not Disturb feature is turned off. If Do Not Disturb is on, you will only see your reminder on the screen.
  - Allow Alarms & reminders permission for the App. This permission is required for Android 12 and above.
1. To add a new reminder, go to the Main Menu and tap **Reminders**. Tap **ADD REMINDER**.
  2. Name your reminder.
  3. Tap the time fields to set the time for the reminder.
  4. Tap **DONE**. You will now see your reminder on the list along with the time you will receive it.

## Note:

- To delete a reminder, swipe the reminder and tap the  symbol. The reminder to check your glucose cannot be deleted.
- Your reminders will be received as notifications that you can swipe or tap to dismiss.

# Settings and Other Options in the Main Menu

## Settings

### App Settings:

**Unit of Measurement** – View the glucose unit of measure used in the App.

**Report Settings** – Work with your healthcare professional to set your Target Glucose Range, which is displayed on glucose graphs in the App and used to calculate your Time In Range. The Target Glucose Range setting will not set glucose alarm levels. Tap **SAVE** when you are done.

**Carbohydrate Units** – Choose grams or portions for food notes that you enter. Tap **SAVE**

when you are done.

**Scan Sounds** – Select whether you would like to hear a sound in addition to a vibration when you scan the Sensor. Remember that Scan Sounds inherit the volume settings on your smartphone. If your smartphone volume is turned off, you will not hear a scan sound. The Scan Sounds setting does not affect alarms. Tap **SAVE** when you are done.

**Text to Speech** – Turn on to have the glucose reading read aloud when you scan the Sensor. You will hear only your current glucose value and trend arrow direction. Additional information, such as the glucose graph and any message, are available on your glucose reading screen. Always review your glucose reading to get complete information. Remember that this feature inherits the volume settings on your phone. If your phone volume is turned off, you will not hear the glucose reading read aloud. Tap **SAVE** when you are done.

**Note:** If you are using a FreeStyle Libre 2 Sensor that you started with the App, you can tap your glucose reading to have it read aloud if this feature is turned on.

## Account Settings:

**Account Details** – View/change your LibreView account information.

**Account Password** – Change your LibreView account password.

## Connected Apps

The **Connected Apps** option in the Main menu opens a web browser within the App. It lists different apps you can connect with to share your data. The available apps may vary based on your country/region. To connect your data with apps listed in this option, select them from the list of apps and follow the onscreen instructions.

## Help

View in-app tutorials, access this User's Manual and review the App's legal information. You can also view a list of events recorded by the App, which may be used by Customer Service to help troubleshoot.

## About

View App software version and other information.

## Activities

**Bathing, Showering and Swimming:** Your Sensor is water-resistant and can be worn while

bathing, showering or swimming. Do NOT take your Sensor deeper than 1 metre (3 ft) or immerse it longer than 30 minutes in water. Note that Bluetooth performance may be impacted if using the System while underwater.

**Sleeping:** Your Sensor should not interfere with your sleep. If you have reminders set to go off while you are sleeping or have glucose alarms set, place your smartphone nearby. The Sensor can store up to 8 hours of data, therefore it is suggested that you review your glucose graph before you go to sleep and when you wake up to check for gaps and scan your Sensor if needed to capture all your data.

**Travelling by Air:** You may use your Sensor while on an aircraft, following any requests from the flight crew.

**IMPORTANT:** Sensor glucose readings and alarms will not be issued while your phone is in aeroplane mode unless Bluetooth is enabled.

- You can continue to scan your Sensor to get glucose readings after you put your phone in aeroplane mode. Confirm NFC is turned on.
- Some airport full-body scanners include X-ray or millimetre radio-wave, which you cannot expose your Sensor to. The effect of these scanners has not been evaluated and the exposure may damage the Sensor or cause inaccurate results. To avoid removing your Sensor, you may request another type of screening. If you do choose to go through a full-body scanner, you must remove your Sensor.
- The Sensor can be exposed to common electrostatic (ESD) and electromagnetic interference (EMI), including airport metal detectors.

**Note:** Changing the time affects the graphs and statistics. The ⌚ symbol may appear on your glucose graph indicating a time change. Gaps in the graph may result or glucose readings may be hidden.

## Maintenance and Disposal

**Maintenance:** The Sensor has no serviceable parts.

**Disposal:** The product contains electronic equipment, batteries, sharps and materials that may contact bodily fluids during use. Dispose of product in accordance with all applicable local regulations. Contact Customer Service for further information on the appropriate disposal of components.

## Troubleshooting

This section lists problems that you may experience, the possible cause(s) and recommended actions. If there is an error, a message will appear on the screen with

directions to resolve the error.

**IMPORTANT:** If you are having issues with the App, please keep in mind that uninstalling the App and/or clearing data will cause you to lose all historical data and end the Sensor currently in use. Please call Customer Service if you have any questions.

## Problems at the Sensor Application Site

Problem: **The Sensor is not sticking to your skin.**

What it may mean: The site is not free of dirt, oil, hair or sweat.

What to do: 1. Remove the Sensor. 2. Clean the site with a plain soap and water and consider shaving. 3. Follow the instructions in [Applying Your Sensor](#) and [Starting Your Sensor](#).

---

Problem: **Skin irritation at the Sensor application site.**

What it may mean: Seams or other constrictive clothing or accessories causing friction at the site **OR** you may be sensitive to the adhesive material.

What to do: Ensure that nothing rubs on the site. If the irritation is where the adhesive touches skin, contact your health care professional to identify the best solution.

## Problems Starting Your Sensor or Receiving Sensor Readings

Display: **Sensor Starting Up**

What it may mean: Sensor is not ready to read glucose.

What to do: Wait until the 60 minute Sensor start-up period has completed.

---

Display: **Signal Loss Alarm**

What it may mean: Sensor has not automatically communicated with the App in the last 20 minutes.

What to do: Make sure your phone is within 6 metres (20 ft) of the Sensor. Try scanning the Sensor to get a glucose reading. Then try turning Bluetooth OFF then ON again. If that doesn't work, try turning your phone OFF then ON again. If **Signal Loss Alarm** shows again, contact Customer Service.

---

Display: **Signal Loss**

What it may mean: Sensor has not automatically communicated with the App in the last 5 minutes.

What to do: Make sure that your phone is within 6 metres (20 ft) of the Sensor and you have not force closed the App. First try scanning your Sensor. Then try turning Bluetooth OFF then

ON again. If that doesn't work, try turning your phone OFF then ON again. If **Signal Loss** shows again, contact Customer Service.

---

Display: **Sensor Ended**

What it may mean: The Sensor life has ended.

What to do: Apply and start a new Sensor.

---

Display: **New Sensor Found**

What it may mean: You scanned a new Sensor before your previous Sensor ended.

What to do: Your phone can only be used with one Sensor at a time. If you start a new Sensor, you will no longer be able to use your old Sensor. If you would like to begin using the new Sensor, select **Yes**.

---

Display: **Sensor Error**

What it may mean: The Sensor is unable to provide a glucose reading.

What to do: Check again after the duration specified in the message.

---

Display: **Glucose Reading Is Unavailable**

What it may mean: The Sensor is unable to provide a glucose reading.

What to do: Check again after the duration specified in the message.

---

Display: **Sensor Too Hot**

What it may mean: Your Sensor is too hot to provide a glucose reading.

What to do: Move to a location where the temperature is appropriate and check again in a few minutes.

---

Display: **Sensor Too Cold**

What it may mean: Your Sensor is too cold to provide a glucose reading.

What to do: Move to a location where the temperature is appropriate and check again in a few minutes.

---

Display: **Check Sensor**

What it may mean: The Sensor tip may not be under your skin.

What to do: Try to start your Sensor again. If you see **Check Sensor** again on the screen,

your Sensor was not applied properly. Apply and start a new Sensor.

---

Display: **Replace Sensor**

What it may mean: The App has detected a problem with your Sensor.

What to do: Apply and start a new Sensor.

---

Display: **Unexpected Application Error**

What it may mean: The App has detected an unexpected error.

What to do: Shut down the App completely and restart.

---

Display: **Incompatible Sensor**

What it may mean: The Sensor cannot be used with the App.

What to do: Call Customer Service.

---

Display: **Scan Error**

What it may mean: Your phone was unable to scan the Sensor **OR** another NFC app is competing for your phone's NFC.

What to do: Try scanning the Sensor again. Make sure to hold your phone still once you get the first tone and/or vibration. Wait until you get the second tone and/or vibration to move it away from the Sensor. Make sure you are not touching any buttons on the smartphone or the screen **OR** open the App to scan your Sensor because another app that uses NFC was detected.

---

Display: **Bluetooth Off**

What it may mean: The Bluetooth setting on your phone is turned off.

What to do: Go to your phone settings and enable Bluetooth.

---

Display: **App Permission Required**

What it may mean: A required App permission is turned off.

What to do: Follow the instructions on the screen to turn the permission on.

## Problems Receiving Glucose Alarms

What it may mean: You have turned glucose alarms off.

What to do: Go to the main menu and then select **Alarms**.

---

What it may mean: The Sensor is not communicating with the App or there may be a problem with the Sensor.

What to do: The Sensor must be within range (6 metres (20 ft)) of your phone for you to receive alarms. Make sure that you are within this range. You will see the  or  symbol when your Sensor has not communicated with the App in 5 minutes. If the **Signal Loss Alarm** is on, you will be notified if there has been no communication for 20 minutes. First try scanning your Sensor. Then try turning Bluetooth OFF then ON again. If that doesn't work, try turning your phone OFF then ON again. If the Signal Loss Alarm persists, contact Customer Service.

---

What it may mean: One or more of the following is turned off: Bluetooth, App Notifications, Nearby Devices permission (Android 12 and above), Alarms & reminders permission (Android 12 and above), Lock Screen notifications, Notification sounds, Channel notifications, Sound and Pop-up notifications, or general phone sounds or vibration. Or, you have turned on Battery Optimisation. Or, you have turned on Do Not Disturb mode without managing the settings for your alarms to override Do Not Disturb.

What to do: Check to make sure that you have the correct settings and permissions enabled on your phone to receive alarms. Go to [Alarms](#) for more information.

---

What it may mean: You may have set an alarm level that is higher or lower than you intended.

What to do: Confirm your alarm settings are appropriate.

---

What it may mean: You have already dismissed this type of alarm.

What to do: You will receive another alarm when a new low or high glucose episode starts.

---

What it may mean: If you are using peripherals such as wireless headphones or a smartwatch you may receive alarms on only one device or peripheral, not all.

What to do: Disconnect headphones or peripherals when you are not using them.

---

What it may mean: FreeStyle Libre 2 app was put to sleep by the phone operating system.

What to do: Put FreeStyle Libre 2 app on the list of apps that will not be put to sleep.

---

What it may mean: You have closed the App.

What to do: Make sure the App is always open in the background.

---

What it may mean: Your Sensor has ended.

What to do: Replace your Sensor with a new one.

## Customer Service

Customer Service is available to answer any questions you may have about FreeStyle Libre 2 app. Please go to [www.FreeStyleLibre.com](http://www.FreeStyleLibre.com) or refer to the product insert in your Sensor Kit for your Customer Service phone number. A printed copy of this User's Manual is available upon request.

## Sensor Specifications

**Sensor glucose assay method:** Amperometric electrochemical sensor

**Sensor glucose reading range:** 2.2 to 22.2 mmol/L

**Sensor size:** 5 mm height and 35 mm diameter

**Sensor weight:** 5 grams

**Sensor power source:** One silver oxide battery

**Sensor life:** Up to 14 days

**Sensor memory:** 8 hours (glucose readings stored every 15 minutes)

**Operating temperature:** 10°C to 45°C

**Sensor Applicator and Sensor Pack storage temperature:** 4°C to 25°C

**Operating and storage relative humidity:** 10-90%, non-condensing

**Sensor water resistance and ingress protection:** IP27: Can withstand immersion into 1 metre (3 ft) of water for up to 30 minutes. Protected against insertion of objects >12 mm diameter.

**Operating and storage altitude:** -381 metres (-1,250 ft) to 3,048 metres (10,000 ft)

**Radio Frequency:** 2.402-2.480 GHz BLE; GFSK; 0dBm EIRP

**Sensor transmission range:** 6 metres (20 ft) unobstructed

# Labelling Symbols and Definitions

Consult instructions for use

---

Temperature limit

---

Manufacturer

---

Batch code

---

Unique Device Identifier

---

Type BF applied part

---

CODE

Sensor code

---

Do not re-use

---

Use-by date

---

Catalogue number

---

Serial number

---

Caution

---

Sterilised by irradiation

---

Humidity limitation

---

Do not use if package is damaged.

---

This product must not be disposed of via municipal waste collection. Contact the manufacturer for details.

---

## Electromagnetic Compatibility

## Model: FreeStyle Libre 2 Sensor

### IC: 12106A-LIB02S

- The Sensor needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in this manual.
- Portable and mobile RF communications equipment can affect the Sensor.
- Use of accessories, transducers and cables other than those specified or provided by Abbott Diabetes Care could result in increased electromagnetic emissions or decreased electromagnetic immunity of the System and result in improper operation.
- The Sensor should not be used adjacent to or stacked with other equipment and that if adjacent or stacked use is necessary, the Sensor should be observed to verify normal operation in the configuration in which it will be used.
- This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.
- Changes or modifications not approved by Abbott could void the user's authority to operate the equipment.

**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Under Industry Canada regulations, this radio transmitter may only operate using an

antenna of a type and maximum (or lesser) gain approved for the transmitter by Industry Canada. To reduce potential radio interference to other users, the antenna type and its gain should be so chosen that the equivalent isotropically radiated power (EIRP) is not more than that necessary for successful communication.

This device complies with Industry Canada's license-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference; and
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

## **Guidance and manufacturer's declaration – electromagnetic emissions**

The Sensor is intended for use in the electromagnetic environment specified below. The customer or the user of the Sensor should ensure that it is used in such an environment.

Emissions test: RF emissions; CISPR 11

Compliance: Group 1

Electromagnetic environment – guidance: The Sensor uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.

Emissions test: RF emissions; CISPR 11

Compliance: Class B

Electromagnetic environment – guidance: The Sensor is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes.

## **Guidance and manufacturer's declaration – electromagnetic immunity**

The Sensor is intended for use in the electromagnetic environment specified below. The customer or the user of the Sensor should assure that it is used in such an environment.

Immunity test: Electrostatic discharge (ESD); IEC 61000-4-2

IEC 60601 test level:  $\pm 8$  kV contact;  $\pm 15$  kV air

Compliance level:  $\pm 8$  kV contact;  $\pm 15$  kV air

Electromagnetic environment – guidance: Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.

Immunity test: Power frequency (50/60 Hz); magnetic field; IEC 61000-4-8

IEC 60601 test level: 30 A/m

Compliance level: 30 A/m

Electromagnetic environment – guidance: Power frequency magnetic fields should be at levels characteristic of a typical location in a typical domestic, commercial or hospital environment.

Immunity test: Radiated RF; IEC 61000-4-3

IEC 60601 test level: 10 V/m; 80 MHz to 2.7 GHz

Compliance level: 10 V/m

Immunity Test: Proximity fields from RF wireless communications equipment; IEC 61000-4-3

Test Level: See the table below

Compliance Level: Compliance to the tested levels

Electromagnetic Environment – guidance:

Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to the Sensor.

Otherwise, degradation of the performance of the System could result.

The table below lists the immunity test levels at specific test frequencies for testing the effects of some wireless communications equipment. The frequencies and services listed in the table are representative examples in health care and in various locations where the System may be used.

| Test frequency (MHz) | Band <sup>a)</sup> (MHz) | Service <sup>a)</sup> | Modulation <sup>b)</sup>                                | Maximum power (W) | Distance (m) | IMMUNITY TEST LEVEL (V/m) |
|----------------------|--------------------------|-----------------------|---|-------------------|--------------|---------------------------|
| 385                  | 380 – 390                | TETRA 400             | Pulse modulation <sup>b)</sup><br>18 Hz                 | 1.8               | 0.3          | 27                        |
| 450                  | 430 – 470                | GMRS 460, FRS 460     | FM <sup>c)</sup><br>$\pm 5$ kHz deviation<br>1 kHz sine | 2                 | 0.3          | 28                        |
|                      |                          |                       |   |                   |              |                           |

|      |             |   |  |     |     |    |
|------|-------------|---|--|-----|-----|----|
| 710  | 704 - 787   | LTE Band 13, 17   | Pulse modulation <sup>b)</sup><br>217 Hz | 0.2 | 0.3 | 9  |
| 745  |             |   |  |     |     |    |
| 780  |             |   |  |     |     |    |
| 810  | 800 - 960   | GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5          | Pulse modulation <sup>b)</sup><br>18 Hz  | 2   | 0.3 | 28 |
| 870  |             |   |  |     |     |    |
| 930  |             |   |  |     |     |    |
| 1720 | 1700 - 1990 | GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS | Pulse modulation <sup>b)</sup><br>217 Hz | 2   | 0.3 | 28 |
| 1845 |             |   |  |     |     |    |
| 1970 |             |   |  |     |     |    |
| 2450 | 2400 - 2570 | Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7            | Pulse modulation <sup>b)</sup><br>217 Hz | 2   | 0.3 | 28 |
| 5240 | 5100 - 5800 | WLAN 802.11 a/n   | Pulse modulation <sup>b)</sup><br>217 Hz | 0.2 | 0.3 | 9  |
| 5500 |             |   |  |     |     |    |
| 5785 |             |   |  |     |     |    |

a) For some services, only the uplink frequencies are included.

b) The carrier shall be modulated using a 50% duty cycle square wave signal.

c) As an alternative to FM modulation, 50% pulse modulation at 18 Hz may be used, because while it does not represent actual modulation, it would be worst case.

Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,<sup>a</sup> should be less than the compliance level in each frequency range.<sup>b</sup>

Interference may occur in the vicinity of equipment marked with the following symbol:

<sup>a</sup> Field strengths from fixed transmitters, such as base stations for radio (mobile/wireless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Sensor is used exceeds the applicable RF compliance level above, the Sensor should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Sensor.

<sup>b</sup> Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 10 V/m.

## Performance Characteristics

### Overview of Clinical Studies

Two studies were conducted in the United States (US) to evaluate the performance, safety, effectiveness and precision of the FreeStyle Libre 2 Flash Glucose Monitoring System (System). One study included adults (Adult study) and the other study included children (Paediatric study). All subjects in both studies required insulin to manage their diabetes. To measure the precision of the System, each subject wore two Sensors, one on the back of each upper arm, for a period of up to 14 days. While in the clinic, subjects had their venous blood glucose analysed using a laboratory reference method, the Yellow Springs Instrument Life Sciences 2300 STAT Plus™ Glucose & Lactate Analyzer (YSI). Sensor glucose readings were then compared to the YSI glucose results to evaluate the System's performance. In the Paediatric study, no YSI measurements were obtained for children under the age of six.

**Adult study:** The Adult study was conducted at 5 centres with 146 subjects in total (91.1% Type 1, 8.9% Type 2), all aged eighteen and older. Subjects had their venous blood glucose analysed over three separate visits to the clinical centre. Each visit lasted up to ten hours. 144 subjects were analysed during the beginning of the Sensor wear period (day 1, 2 or 3), 91 subjects were analysed during the early middle period (day 7 or 8), 55 subjects were analysed during the late middle period (day 9 or 12) and 76 subjects were analysed during the end period (day 13 or 14). During each visit, adult subjects had their glucose levels deliberately manipulated per the study protocol to raise or lower glucose. This was done to assess performance of the System over the range that the System measures glucose (2.2 – 22.2 mmol/L).

**Paediatric study:** The Paediatric study was conducted at 4 centres with 139 subjects in total (98.6% Type 1, 1.4% Type 2), all aged four to seventeen. Subjects age six and older had their venous blood glucose analysed for up to 16 hours over one or two separate visits to the clinical centre. Each visit lasted up to eight hours. During each visit, subjects age 11 and older had their glucose levels deliberately manipulated per the study protocol to raise or lower glucose. This was done to assess performance of the System over the range that the System measures glucose (2.2 – 22.2 mmol/L). 48 subjects were analysed during the beginning of the Sensor wear period (day 1 or 2), 50 subjects were analysed during the early middle period (day 7 or 8), 51 subjects were analysed during the late middle period (day 9 or 12) and 51 subjects were analysed during the end period (day 13 or 14). All subjects tested their blood glucose using fingerstick capillary samples at least four times during each day of the study.

## **Accuracy**

Accuracy of the System was measured by comparing paired System Glucose Measurement (GM) and YSI blood glucose values. The percentage of total System readings that were within  $\pm 1.11$  mmol/L for YSI blood glucose values  $< 3.9$  mmol/L or  $\pm 20\%$  of YSI for blood glucose values  $\geq 3.9$  mmol/L is displayed in [Table 1](#). The Mean Absolute Relative Difference (MARD) gives an indication of the average percent disagreement between the GM and the reference. For example, in the Adult study, 92.3% of the readings fell within  $\pm 1.11$  mmol/L of YSI blood

| Glucose Level† (mmol/L) | of GM-Reference Pairs | Percent Within $\pm 0.83$ mmol/L | Percent Within $\pm 1.11$ mmol/L | Percent Within $\pm 2.2$ mmol/L | Percent Within $\pm 15\%$ | Percent Within $\pm 20\%$ | Percent Within $\pm 40\%$ | Mean bias (mmol/L) | MARD (%) |
|-------------------------|-----------------------|----------------------------------|----------------------------------|---------------------------------|---------------------------|---------------------------|---------------------------|--------------------|----------|
| <3.0                    | 139                   | 71.9                             | 79.1                             | 97.1                            |                           |                           |                           | -0.55              | 17.1     |
| 3.0 – 3.9               | 863                   | 86.4                             | 90.5                             | 97.1                            |                           |                           |                           | -0.27              | 12.0     |
| 3.9 – 10.0              | 2690                  |                                  |                                  |                                 | 77.4                      | 87.6                      | 98.7                      | -0.19              | 10.6     |
| 10.0 – 13.9             | 1236                  |                                  |                                  |                                 | 86.0                      | 94.7                      | 99.7                      | -0.50              | 8.3      |
| 13.9 – 19.4             | 1411                  |                                  |                                  |                                 | 91.6                      | 97.6                      | 99.9                      | -0.23              | 7.3      |
| >19.4                   | 207                   |                                  |                                  |                                 | 96.6                      | 98.6                      | 98.6                      | 0.63               | 6.6      |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

† System range is 2.2 – 22.2 mmol/L.

**Table 4.** Accuracy to YSI within YSI Glucose Ranges (Adult; n=144)

| YSI Glucose Level (mmol/L) | Number of GM-Reference Pairs | Percent Within $\pm 0.83$ mmol/L | Percent Within $\pm 1.11$ mmol/L | Percent Within $\pm 2.2$ mmol/L | Percent Within $\pm 15\%$ | Percent Within $\pm 20\%$ | Percent Within $\pm 40\%$ | Mean bias (mmol/L) | MARD (%) |
|----------------------------|------------------------------|----------------------------------|----------------------------------|---------------------------------|---------------------------|---------------------------|---------------------------|--------------------|----------|
| <3.0                       | 440                          | 91.1                             | 97.5                             | 100.0                           |                           |                           |                           | 0.41               | 15.5     |
| 3.0 – 3.9                  | 3033                         | 94.7                             | 98.6                             | 100.0                           |                           |                           |                           | 0.08               | 10.2     |
| 3.9 – 10.0                 | 7544                         |                                  |                                  |                                 | 77.3                      | 86.8                      | 99.4                      | -0.27              | 10.4     |
| 10.0 – 13.9                | 2950                         |                                  |                                  |                                 | 87.7                      | 93.5                      | 99.6                      | -0.45              | 8.1      |
| 13.9 – 19.4                | 4142                         |                                  |                                  |                                 | 91.8                      | 96.3                      | 99.6                      | -0.54              | 6.8      |
| >19.4                      | 698                          |                                  |                                  |                                 | 85.2                      | 93.4                      | 100.0                     | -1.32              | 7.9      |

**Table 5.** Accuracy to YSI within YSI Glucose Ranges (Paediatric\*; n=129)

| YSI Glucose Level (mmol/L) | Number of GM-Referenced Pairs | Percent Within $\pm 0.83$ mmol/L | Percent Within $\pm 1.11$ mmol/L | Percent Within $\pm 2.2$ mmol/L | Percent Within $\pm 15\%$ | Percent Within $\pm 20\%$ | Percent Within $\pm 40\%$ | Mean bias (mmol/L) | MARD (%) |
|----------------------------|-------------------------------|----------------------------------|----------------------------------|---------------------------------|---------------------------|---------------------------|---------------------------|--------------------|----------|
| <3.0                       | 131                           | 93.9                             | 98.5                             | 100.0                           |                           |                           |                           | 0.37               | 14.2     |
| 3.0 – 3.9                  | 751                           | 96.5                             | 98.8                             | 100.0                           |                           |                           |                           | 0.06               | 9.3      |
| 3.9 – 10.0                 | 2743                          |                                  |                                  |                                 | 74.3                      | 84.8                      | 98.0                      | -0.16              | 11.4     |
| 10.0 – 13.9                | 1104                          |                                  |                                  |                                 | 86.6                      | 92.9                      | 99.0                      | -0.21              | 8.4      |
| 13.9 – 19.4                | 1572                          |                                  |                                  |                                 | 90.1                      | 97.5                      | 99.9                      | -0.49              | 7.7      |
| >19.4                      | 245                           |                                  |                                  |                                 | 90.6                      | 97.6                      | 100.0                     | -1.06              | 7.1      |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

## Agreement with 'LO' and 'HI' GM Reading against YSI Reference

The System reports glucose concentrations between 2.2 and 22.2 mmol/L. When the System determines that glucose level is below 2.2 mmol/L, it will report as 'LO'. When the System determines that glucose level is above 22.2 mmol/L, it will report as 'HI'. [Table 6](#) and [Table 7](#) display the concurrence between the GM and YSI reference glucose when GM reads 'LO'. For example, in the Adult study, when GM reading was 'LO', YSI glucose values were less than 2.8 mmol/L 20.0% of the time, less than 3.3 mmol/L 40.0% of the time, less than 3.9 mmol/L 40.0% of the time, less than 4.4 mmol/L 80.0% of the time and equal to or above 4.4 mmol/L 20.0% of the time.

**Table 6.** Concurrence Analysis with 'LO' GM Reading (Adult; n=144)

| GM - Reference Pairs | YSI (mmol/L) |      |      |      |      | N |
|----------------------|--------------|------|------|------|------|---|
|                      | <2.8         | <3.3 | <3.9 | <4.4 | ≥4.4 |   |
| n                    | 1            | 2    | 2    | 4    | 1    | 5 |
| Cumulative %         | 20.0         | 40.0 | 40.0 | 80.0 | 20.0 |   |

**Table 7.** Concurrence Analysis with 'LO' GM Reading (Paediatric\*; n=129)

| GM - Reference Pairs | YSI (mmol/L) |      |       |       |      | N |
|----------------------|--------------|------|-------|-------|------|---|
|                      | <2.8         | <3.3 | <3.9  | <4.4  | ≥4.4 |   |
| n                    | 0            | 1    | 2     | 2     | 0    | 2 |
| Cumulative %         | 0.0          | 50.0 | 100.0 | 100.0 | 0.0  |   |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

[Table 8](#) and [Table 9](#) display the concurrence between the GM and YSI reference glucose when GM reads 'HI'. In the Adult study, when GM reading was 'HI', YSI glucose values were above 22.2 mmol/L 34.4% of the time, above 19.4 mmol/L 97.5% of the time, above 16.7 mmol/L 100.0% of the time, above 13.9 mmol/L 100.0% of the time and below or equal to 13.9 mmol/L 0.0% of the time.

**Table 8.** Concurrence Analysis with 'HI' GM Reading (Adult; n=144)

| GM-Reference Pairs | YSI (mmol/L) |       |       |       |       | N   |
|--------------------|--------------|-------|-------|-------|-------|-----|
|                    | >22.2        | >19.4 | >16.7 | >13.9 | ≤13.9 |     |
| n                  | 42           | 119   | 122   | 122   | 0     | 122 |
| Cumulative %       | 34.4         | 97.5  | 100.0 | 100.0 | 0.0   |     |

**Table 9.** Concurrence Analysis with 'HI' GM Reading (Paediatric\*; n=129)

| GM-Reference Pairs | YSI (mmol/L) |       |       |       |       | N  |
|--------------------|--------------|-------|-------|-------|-------|----|
|                    | >22.2        | >19.4 | >16.7 | >13.9 | ≤13.9 |    |
| n                  | 25           | 40    | 43    | 45    | 0     | 45 |
| Cumulative %       | 55.6         | 88.9  | 95.6  | 100.0 | 0.0   |    |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

## Concurrence of System and Reference (GM vs. YSI)

The percentage of concurring glucose values (GM vs. YSI) in each glucose reference range is presented for each GM range in [Table 10](#) and [Table 11](#) and for each YSI range in [Table 12](#) and [Table 13](#). For example, in the Adult study, when the System glucose readings were within the 4.5 and 6.7 mmol/L range, actual blood glucose values were between 2.2 and 3.3 mmol/L 0.2% of the time, between 3.4 and 4.4 mmol/L 10.9% of the time, between 4.5 and 6.7 mmol/L 70.0% of the time, between 6.7 and 8.9 mmol/L 17.9% of the time, between 8.9 and 11.1 mmol/L 0.9% of the time, and between 11.2 and 13.9 mmol/L 0.1% of the time.

**Table 10.** Concurrence Analysis by GM Glucose Level (Adult; n=144)

| GM Glucose Level (mmol/L) | YSI Glucose Level (mmol/L) |           |           |           |           |            |             |             |             |             |       | N    |
|---------------------------|----------------------------|-----------|-----------|-----------|-----------|------------|-------------|-------------|-------------|-------------|-------|------|
|                           | <2.2                       | 2.2 – 3.3 | 3.4 – 4.4 | 4.5 – 6.7 | 6.7 – 8.9 | 8.9 – 11.1 | 11.2 – 13.9 | 13.9 – 16.7 | 16.7 – 19.4 | 19.5 – 22.2 | >22.2 |      |
| <2.2 <sup>†</sup>         | 20.0                       | 20.0      | 40.0      | 20.0      | -         | -          | -           | -           | -           | -           | -     | 5    |
| 2.2 – 3.3                 | 0.4                        | 52.7      | 43.6      | 3.3       | -         | 0.1        | -           | -           | -           | -           | -     | 1899 |
| 3.4 – 4.4                 | -                          | 18.8      | 62.6      | 18.2      | 0.4       | 0.0        | -           | -           | -           | -           | -     | 3095 |
| 4.5 – 6.7                 | -                          | 0.2       | 10.9      | 70.0      | 17.9      | 0.9        | 0.1         | -           | -           | -           | -     | 3055 |
| 6.7 – 8.9                 | -                          | -         | 0.1       | 9.1       | 69.6      | 19.1       | 1.6         | 0.3         | 0.2         | -           | -     | 2425 |
| 8.9 – 11.1                | -                          | -         | -         | -         | 10.6      | 60.5       | 27.0        | 1.6         | 0.3         | -           | -     | 1753 |
| 11.2 – 13.9               | -                          | -         | -         | -         | -         | 7.0        | 65.4        | 25.6        | 1.9         | 0.1         | -     | 2184 |
| 13.9 – 16.7               | -                          | -         | -         | -         | -         | 0.1        | 8.4         | 66.8        | 22.8        | 1.9         | 0.1   | 2330 |
| 16.7 – 19.4               | -                          | -         | -         | -         | -         | -          | 0.4         | 13.6        | 68.8        | 16.0        | 1.2   | 1528 |
| 19.5 – 22.2               | -                          | -         | -         | -         | -         | -          | -           | 0.6         | 27.3        | 63.6        | 8.6   | 538  |
| >22.2 <sup>†</sup>        | -                          | -         | -         | -         | -         | -          | -           | -           | 2.5         | 63.1        | 34.4  | 122  |

<sup>†</sup> Levels out of System dynamic range.

**Table 11.** Concurrence Analysis by GM Glucose Level (Paediatric\*; n=129)

| GM Glucose Level (mmol/L) | YSI Glucose Level (mmol/L) |           |           |           |           |            |             |             |             |             |       | N    |
|---------------------------|----------------------------|-----------|-----------|-----------|-----------|------------|-------------|-------------|-------------|-------------|-------|------|
|                           | <2.2                       | 2.2 – 3.3 | 3.4 – 4.4 | 4.5 – 6.7 | 6.7 – 8.9 | 8.9 – 11.1 | 11.2 – 13.9 | 13.9 – 16.7 | 16.7 – 19.4 | 19.5 – 22.2 | >22.2 |      |
| <2.2 <sup>†</sup>         | -                          | 50.0      | 50.0      | -         | -         | -          | -           | -           | -           | -           | -     | 2    |
| 2.2 – 3.3                 | 0.6                        | 48.6      | 42.5      | 7.8       | 0.6       | -          | -           | -           | -           | -           | -     | 527  |
| 3.4 – 4.4                 | -                          | 12.1      | 61.9      | 24.3      | 1.7       | -          | -           | -           | -           | -           | -     | 915  |
| 4.5 – 6.7                 | -                          | 0.2       | 11.2      | 69.0      | 18.2      | 1.3        | 0.1         | -           | -           | -           | -     | 1006 |
| 6.7 – 8.9                 | -                          | -         | -         | 11.4      | 71.0      | 15.8       | 1.8         | -           | -           | -           | -     | 868  |

|                    |   |   |   |     |      |      |      |      |      |      |      |     |
|--------------------|---|---|---|-----|------|------|------|------|------|------|------|-----|
| 8.9 – 11.1         | - | - | - | 0.1 | 18.2 | 61.3 | 20.1 | 0.3  | -    | -    | -    | 703 |
| 11.2 – 13.9        | - | - | - | -   | 0.2  | 9.6  | 55.3 | 33.6 | 1.2  | 0.1  | -    | 909 |
| 13.9 – 16.7        | - | - | - | -   | -    | 0.1  | 14.1 | 60.8 | 23.7 | 1.3  | -    | 818 |
| 16.7 – 19.4        | - | - | - | -   | -    | -    | 0.3  | 24.8 | 58.2 | 16.5 | 0.2  | 593 |
| 19.5 – 22.2        | - | - | - | -   | -    | 1.0  | -    | 0.5  | 33.8 | 59.4 | 5.3  | 207 |
| >22.2 <sup>†</sup> | - | - | - | -   | -    | -    | -    | 4.4  | 6.7  | 33.3 | 55.6 | 45  |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

<sup>†</sup> Levels out of System dynamic range.

**Table 12.** Concurrence Analysis by YSI Glucose Level (Adult; n=144)

| YSI Glucose Level (mmol/L) | GM Glucose Level (mmol/L) |           |           |           |           |            |             |             |             |             |                    | N    |
|----------------------------|---------------------------|-----------|-----------|-----------|-----------|------------|-------------|-------------|-------------|-------------|--------------------|------|
|                            | <2.2 <sup>†</sup>         | 2.2 – 3.3 | 3.4 – 4.4 | 4.5 – 6.7 | 6.7 – 8.9 | 8.9 – 11.1 | 11.2 – 13.9 | 13.9 – 16.7 | 16.7 – 19.4 | 19.5 – 22.2 | >22.2 <sup>†</sup> |      |
| <2.2                       | 12.5                      | 87.5      | -         | -         | -         | -          | -           | -           | -           | -           | -                  | 8    |
| 2.2 – 3.3                  | 0.1                       | 62.9      | 36.6      | 0.4       | -         | -          | -           | -           | -           | -           | -                  | 1591 |
| 3.4 – 4.4                  | 0.1                       | 26.7      | 62.4      | 10.8      | 0.1       | -          | -           | -           | -           | -           | -                  | 3103 |
| 4.5 – 6.7                  | 0.0                       | 2.1       | 18.9      | 71.6      | 7.4       | -          | -           | -           | -           | -           | -                  | 2984 |
| 6.7 – 8.9                  | -                         | -         | 0.5       | 22.5      | 69.5      | 7.6        | -           | -           | -           | -           | -                  | 2430 |
| 8.9 – 11.1                 | -                         | 0.1       | 0.1       | 1.6       | 27.1      | 62.1       | 9.0         | 0.1         | -           | -           | -                  | 1709 |
| 11.2 – 13.9                | -                         | -         | -         | 0.1       | 1.9       | 22.1       | 66.6        | 9.1         | 0.3         | -           | -                  | 2147 |
| 13.9 – 16.7                | -                         | -         | -         | -         | 0.3       | 1.2        | 23.7        | 65.9        | 8.8         | 0.1         | -                  | 2361 |
| 16.7 – 19.4                | -                         | -         | -         | -         | 0.3       | 0.3        | 2.3         | 29.8        | 59.0        | 8.2         | 0.2                | 1784 |
| 19.5 – 22.2                | -                         | -         | -         | -         | -         | -          | 0.3         | 6.2         | 34.4        | 48.2        | 10.9               | 709  |
| >22.2                      | -                         | -         | -         | -         | -         | -          | -           | 1.9         | 16.7        | 42.6        | 38.9               | 108  |

<sup>†</sup> Levels out of System dynamic range.

**Table 13.** Concurrence Analysis by YSI Glucose Level (Paediatric\*; n=129)

| YSI Glucose Level (mmol/L) | GM Glucose Level (mmol/L) |           |           |           |           |            |             |             |             |             |                    | N    |
|----------------------------|---------------------------|-----------|-----------|-----------|-----------|------------|-------------|-------------|-------------|-------------|--------------------|------|
|                            | <2.2 <sup>†</sup>         | 2.2 – 3.3 | 3.4 – 4.4 | 4.5 – 6.7 | 6.7 – 8.9 | 8.9 – 11.1 | 11.2 – 13.9 | 13.9 – 16.7 | 16.7 – 19.4 | 19.5 – 22.2 | >22.2 <sup>†</sup> |      |
| <2.2                       | -                         | 100.0     | -         | -         | -         | -          | -           | -           | -           | -           | -                  | 3    |
| 2.2 – 3.3                  | 0.3                       | 69.2      | 30.0      | 0.5       | -         | -          | -           | -           | -           | -           | -                  | 370  |
| 3.4 – 4.4                  | 0.1                       | 24.8      | 62.6      | 12.5      | -         | -          | -           | -           | -           | -           | -                  | 904  |
| 4.5 – 6.7                  | -                         | 3.9       | 21.0      | 65.7      | 9.4       | 0.1        | -           | -           | -           | -           | -                  | 1057 |
| 6.7 – 8.9                  | -                         | 0.3       | 1.7       | 19.3      | 65.0      | 13.5       | 0.2         | -           | -           | -           | -                  | 948  |

|             |   |   |   |     |      |      |      |      |      |      |      |     |
|-------------|---|---|---|-----|------|------|------|------|------|------|------|-----|
| 8.9 – 11.1  | - | - | - | 1.9 | 20.4 | 64.2 | 13.0 | 0.1  | -    | 0.3  | -    | 671 |
| 11.2 – 13.9 | - | - | - | 0.1 | 2.1  | 18.1 | 64.7 | 14.8 | 0.3  | -    | -    | 778 |
| 13.9 – 16.7 | - | - | - | -   | -    | 0.2  | 32.0 | 52.1 | 15.4 | 0.1  | 0.2  | 954 |
| 16.7 – 19.4 | - | - | - | -   | -    | -    | 1.8  | 31.1 | 55.4 | 11.2 | 0.5  | 623 |
| 19.5 – 22.2 | - | - | - | -   | -    | -    | 0.4  | 4.4  | 39.5 | 49.6 | 6.0  | 248 |
| >22.2       | - | - | - | -   | -    | -    | -    | -    | 2.7  | 29.7 | 67.6 | 37  |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

† Levels out of System dynamic range.

### Glucose Rate of Change Accuracy

The System's glucose rate of change (ROC) accuracy, as assessed by concurrence analysis, is presented in [Table 14](#) and [Table 15](#). For example, in the Adult study, when the Sensor glucose rate of change indicated that glucose was changing slowly (-0.06 to 0.06 mmol/L/min), actual glucose levels in the body were falling quickly (<-0.1 mmol/L/min) 1.1% of the time, falling (-0.1 to -0.06 mmol/L/min) 6.5% of the time, changing slowly (-0.06 to 0.06 mmol/L/min) 81.7% of the time, rising (0.06 to 0.1 mmol/L/min) 7.3% of the time and were rising quickly (>0.1 mmol/L/min) 3.5% of the time.

**Table 14.** Concurrence Analysis by Glucose Rate of Change (Adult; n=144)

| GM Rate (mmol/L/min) | YSI Rate (mmol/L/min) |               |               |             |      | N     |
|----------------------|-----------------------|---------------|---------------|-------------|------|-------|
|                      | <-0.1                 | [-0.1, -0.06] | [-0.06, 0.06] | (0.06, 0.1] | >0.1 |       |
| <-0.1                | 33.8                  | 44.5          | 21.3          | 0.3         | -    | 328   |
| [-0.1 to -0.06]      | 6.8                   | 46.5          | 45.2          | 0.9         | 0.5  | 1092  |
| [-0.06 to 0.06]      | 1.1                   | 6.5           | 81.7          | 7.3         | 3.5  | 14865 |
| (0.06 to 0.1]        | 0.1                   | 1.6           | 39.3          | 38.4        | 20.5 | 1154  |
| >0.1                 | 0.1                   | 0.2           | 17.7          | 32.9        | 49.0 | 885   |

**Table 15.** Concurrence Analysis by Glucose Rate of Change (Paediatric\*; n=129)

| GM Rate (mmol/L/min) | YSI Rate (mmol/L/min) |               |               |             |      | N    |
|----------------------|-----------------------|---------------|---------------|-------------|------|------|
|                      | <-0.1                 | [-0.1, -0.06] | [-0.06, 0.06] | (0.06, 0.1] | >0.1 |      |
| <-0.1                | 44.1                  | 44.7          | 11.2          | -           | -    | 170  |
| [-0.1 to -0.06]      | 11.4                  | 49.5          | 38.0          | 0.4         | 0.6  | 463  |
| [-0.06 to 0.06]      | 1.8                   | 8.7           | 75.2          | 8.8         | 5.4  | 4682 |
| (0.06 to 0.1]        | 0.2                   | 2.6           | 40.2          | 35.5        | 21.5 | 498  |
| >0.1                 | -                     | 0.9           | 19.2          | 29.7        | 50.2 | 448  |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

## Alarm Performance

The tables in this section show the accuracy of the System's Low and High Glucose Alarms. The Alarm Rate tells you how often the alarm is right or wrong. The Detection Rate tells you how often the System is able to recognise and notify you about a low or high glucose event.

### Low Glucose Alarm Performance

[Table 16](#) and [Table 17](#) display the percentages for these parameters:

#### *True Alarm Rate*

Tells you: When you got a low glucose alarm, were you actually low?

Definition: Percentage of time the alarm issued and blood glucose was below the alarm level within 15 minutes before or after the alarm.

#### *False Alarm Rate*

Tells you: Did you get a low glucose alarm that you shouldn't have?

Definition: Percentage of time the alarm issued and blood glucose was not below the alarm level within 15 minutes before or after the alarm.

#### *Detection Rate*

Tells you: When you were low, did you get a low glucose alarm?

Definition: Percentage of time blood glucose was below the alarm level and the alarm issued within 15 minutes before or after the glucose event.

#### *Missed Detection Rate*

Tells you: When you were low, did you miss a low glucose alarm?

Definition: Percentage of time blood glucose was below the alarm level and the alarm didn't issue within 15 minutes before or after the glucose event.

For example, the Adult study found that for a Low Glucose alarm level set to 3.9 mmol/L:

86.0% of the time a low glucose alarm was received when blood glucose was indeed below the alarm level but 14.0% of the time a low glucose alarm was received when blood glucose wasn't below the alarm level.

89.3% of the time blood glucose was below the alarm level and a low glucose alarm was appropriately issued but 10.7% of the time the glucose event was missed and no alarm was issued.

**Table 16.** Low Glucose Alarm Performance (Adult; n=144)

| Low Glucose Alarm Level (mmol/L) | Alarm Rate          |                      | Detection Rate             |                           |
|----------------------------------|---------------------|----------------------|----------------------------|---------------------------|
|                                  | True Alarm Rate (%) | False Alarm Rate (%) | Correct Detection Rate (%) | Missed Detection Rate (%) |
| 3.3                              | 72.3                | 27.7                 | 75.7                       | 24.3                      |
| 3.9                              | 86.0                | 14.0                 | 89.3                       | 10.7                      |
| 4.4                              | 91.3                | 8.7                  | 97.3                       | 2.7                       |
| 5.0                              | 93.6                | 6.4                  | 98.5                       | 1.5                       |

**Table 17.** Low Glucose Alarm Performance (Paediatric<sup>\*</sup>; n=129)

| Low Glucose Alarm Level (mmol/L) | Alarm Rate          |                      | Detection Rate             |                           |
|----------------------------------|---------------------|----------------------|----------------------------|---------------------------|
|                                  | True Alarm Rate (%) | False Alarm Rate (%) | Correct Detection Rate (%) | Missed Detection Rate (%) |
| 3.3                              | 62.9                | 37.1                 | 87.4                       | 12.6                      |
| 3.9                              | 80.3                | 19.7                 | 93.5                       | 6.5                       |
| 4.4                              | 85.6                | 14.4                 | 96.4                       | 3.6                       |
| 5.0                              | 92.2                | 7.8                  | 97.3                       | 2.7                       |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

## High Glucose Alarm Performance

Table 18 and Table 19 display the percentages for these parameters:

### *True Alarm Rate*

Tells you: When you got a high glucose alarm, were you actually high?

Definition: Percentage of time the alarm issued and blood glucose was above the alarm level within 15 minutes before or after the alarm.

### *False Alarm Rate*

Tells you: Did you get a high glucose alarm that you shouldn't have?

Definition: Percentage of time the alarm issued and blood glucose was not above the alarm level within 15 minutes before or after the alarm.

### *Detection Rate*

Tells you: When you were high, did you get a high glucose alarm?

Definition: Percentage of time blood glucose was above the alarm level and the alarm issued within 15 minutes before or after the glucose event.

### Missed Detection Rate

Tells you: When you were high, did you miss a high glucose alarm?

Definition: Amount of time blood glucose was above the alarm level and the alarm didn't issue within 15 minutes before or after the glucose event.

For example, the Adult study found that for a High Glucose alarm level set to 11.1 mmol/L:

99.2% of the time a high glucose alarm was received when blood glucose was indeed above the alarm level but 0.8% of the time a high glucose alarm was received when blood glucose wasn't actually above the alarm level.

97.0% of the time blood glucose was above the alarm level and a high glucose alarm was appropriately issued but 3.0% of the time the glucose event was missed and no alarm was issued.

**Table 18.** High Glucose Alarm Performance (Adult; n=144)

| High Glucose Alarm level (mmol/L) | Alarm Rate          |                      | Detection Rate             |                           |
|-----------------------------------|---------------------|----------------------|----------------------------|---------------------------|
|                                   | True Alarm Rate (%) | False Alarm Rate (%) | Correct Detection Rate (%) | Missed Detection Rate (%) |
| 6.7                               | 99.1                | 0.9                  | 98.2                       | 1.8                       |
| 7.8                               | 99.1                | 0.9                  | 98.1                       | 1.9                       |
| 10.0                              | 99.2                | 0.8                  | 97.8                       | 2.2                       |
| 11.1                              | 99.2                | 0.8                  | 97.0                       | 3.0                       |
| 12.2                              | 99.0                | 1.0                  | 96.9                       | 3.1                       |
| 13.3                              | 98.4                | 1.6                  | 95.6                       | 4.4                       |
| 16.7                              | 96.3                | 3.7                  | 90.0                       | 10.0                      |

**Table 19.** High Glucose Alarm Performance (Paediatric\*; n=129)

| High Glucose Alarm level (mmol/L) | Alarm Rate          |                      | Detection Rate             |                           |
|-----------------------------------|---------------------|----------------------|----------------------------|---------------------------|
|                                   | True Alarm Rate (%) | False Alarm Rate (%) | Correct Detection Rate (%) | Missed Detection Rate (%) |
| 6.7                               | 98.8                | 1.2                  | 98.2                       | 1.8                       |
| 7.8                               | 98.0                | 2.0                  | 98.4                       | 1.6                       |
| 10.0                              | 98.4                | 1.6                  | 98.0                       | 2.0                       |
| 11.1                              | 98.0                | 2.0                  | 98.0                       | 2.0                       |
| 12.2                              | 98.2                | 1.8                  | 96.9                       | 3.1                       |
| 13.3                              | 98.0                | 2.0                  | 95.7                       | 4.3                       |
| 16.7                              | 90.8                | 9.2                  | 91.0                       | 9.0                       |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

## Sensor Accuracy Over Time

The Sensor can be worn for up to 14 days. The percentage of System readings within YSI values and the Mean Absolute Relative Difference (MARD) is presented for the following different wear periods in [Table 20](#) and [Table 21](#): Beginning (Adult: 144 Subjects, Day 1, 2 or 3; Paediatric: 48 Subjects, Day 1 or 2) Early Middle (Adult: 91 Subjects, Day 7 or 8; Paediatric: 50 Subjects, Day 7 or 8), Late Middle (Adult: 55 Subjects, Day 9 or 12; Paediatric: 51 Subjects, Day 9 or 12), and End (Adult: 76 Subjects, Day 13 or 14; Paediatric: 51 Subjects, Day 13 or 14). For values  $\geq 3.9$  mmol/L, the percentage of readings within 15%, 20% and 40% of the YSI value was calculated. For values  $< 3.9$  mmol/L, the percentage of readings within 0.83 mmol/L, 1.11 mmol/L and 2.2 mmol/L of the YSI value was calculated. The percentage of System readings within YSI values for the first 12 hours of a Sensor wear is presented in [Table 22](#).

**Table 20.** Sensor Accuracy Relative to YSI Over the Wear Duration (Adult; n=144)

| Wear Period  | Number of GM-reference pairs | MARD (%) | Within $\pm 15\%$ / $\pm 0.83$ mmol/L | Within $\pm 20\%$ / $\pm 1.11$ mmol/L | Within $\pm 40\%$ / $\pm 2.2$ mmol/L |
|--------------|------------------------------|----------|---------------------------------------|---------------------------------------|--------------------------------------|
| Beginning    | 7027                         | 10.0     | 83.0                                  | 90.2                                  | 99.3                                 |
| Early Middle | 4522                         | 8.5      | 87.7                                  | 94.5                                  | 99.8                                 |
| Late Middle  | 3503                         | 8.8      | 86.8                                  | 93.4                                  | 99.7                                 |
| End          | 3755                         | 9.1      | 86.4                                  | 92.9                                  | 100.0                                |

**Table 21.** Sensor Accuracy Relative to YSI Over the Wear Duration (Paediatric<sup>\*</sup>; n=129)

| Wear Period  | Number of GM-reference pairs | MARD (%) | Within $\pm 15\%$ / $\pm 0.83$ mmol/L | Within $\pm 20\%$ / $\pm 1.11$ mmol/L | Within $\pm 40\%$ / $\pm 2.2$ mmol/L |
|--------------|------------------------------|----------|---------------------------------------|---------------------------------------|--------------------------------------|
| Beginning    | 1828                         | 10.7     | 79.6                                  | 88.5                                  | 98.6                                 |
| Early Middle | 1642                         | 8.0      | 89.5                                  | 94.2                                  | 98.5                                 |
| Late Middle  | 1534                         | 9.7      | 83.6                                  | 92.9                                  | 99.5                                 |
| End          | 1542                         | 10.2     | 82.6                                  | 91.1                                  | 99.3                                 |

\* Includes children 6 – 17 years of age. No YSI measurements were obtained for children 4 – 5 years of age.

**Table 22.** Sensor Accuracy Relative to YSI over the First 12 Hours of Sensor Wear

| Population | Number of GM-reference pairs | Percent Within $\pm 20\%$ / $\pm 1.11$ mmol/L | Percent Within $\pm 40\%$ / $\pm 2.2$ mmol/L |
|------------|------------------------------|---|--|
| Adult      | 1292                         | 81.0  | 97.4   |
| Paediatric | 548                          | 80.3  | 97.3   |

## Sensor Wear Duration

The Sensor can be worn for up to 14 days. To estimate how long a Sensor will work over the wear duration, 146 Sensors were evaluated in the Adult study and 139 Sensors were

evaluated in the Paediatric study to determine how many days of readings each Sensor provided. Of the 146 Sensors in the Adult study, 71.1% lasted until the final day of use. In the Paediatric study, 78.1% of the Sensors lasted until the final day of use. [Table 23](#) and [Table 24](#) display the data for each day in the wear duration for the Adult and Paediatric studies.

**Table 23.** Sensor Survival Rate Over Wear Duration (Adult; n=146)

| Day of Wear | Number of Sensors | Survival Rate (%) |
|-------------|-------------------|-------------------|
| 1           | 145               | 99.3              |
| 2           | 142               | 97.3              |
| 3           | 140               | 95.9              |
| 4           | 137               | 93.8              |
| 5           | 134               | 91.8              |
| 6           | 133               | 91.1              |
| 7           | 132               | 90.4              |
| 8           | 127               | 87.0              |
| 9           | 123               | 84.9              |
| 10          | 119               | 82.2              |
| 11          | 112               | 77.3              |
| 12          | 111               | 76.6              |
| 13          | 104               | 71.8              |
| 14          | 100               | 71.1              |

**Table 24.** Sensor Survival Rate Over Wear Duration (Paediatric; n=139)

| Day of Wear | Number of Sensors | Survival Rate (%) |
|-------------|-------------------|-------------------|
| 1           | 137               | 98.6              |
| 2           | 136               | 97.8              |
| 3           | 134               | 97.1              |
| 4           | 133               | 96.4              |
| 5           | 133               | 96.4              |
| 6           | 133               | 96.4              |
| 7           | 133               | 96.4              |
| 8           | 131               | 94.9              |
| 9           | 126               | 91.3              |
| 10          | 124               | 89.9              |
| 11          | 122               | 88.4              |
| 12          | 120               | 87.0              |
| 13          | 114               | 83.4              |

|    |     |      |
|----|-----|------|
| 14 | 104 | 78.1 |
|----|-----|------|

## Glucose Reading Availability

The System is designed to show a Sensor glucose reading after each scan that is performed throughout the wear period after the start-up time. [Table 25](#) and [Table 26](#) show the glucose reading capture rate for each day of the wear duration.

**Table 25.** Glucose Reading Capture Rate Over Wear Duration (Adult; n=146)

| Day of Wear | Number of Sensors | Capture Rate (%) |
|-------------|-------------------|------------------|
| 1           | 146               | 98.2             |
| 2           | 145               | 98.0             |
| 3           | 143               | 98.2             |
| 4           | 140               | 98.3             |
| 5           | 138               | 98.3             |
| 6           | 135               | 98.3             |
| 7           | 134               | 98.3             |
| 8           | 131               | 98.4             |
| 9           | 128               | 98.4             |
| 10          | 123               | 98.4             |
| 11          | 120               | 98.4             |
| 12          | 113               | 98.5             |
| 13          | 112               | 98.5             |
| 14          | 104               | 98.5             |

**Table 26.** Glucose Reading Capture Rate Over Wear Duration (Paediatric; n=139)

| Day of Wear | Number of Sensors | Capture Rate (%) |
|-------------|-------------------|------------------|
| 1           | 139               | 94.4             |
| 2           | 137               | 94.7             |
| 3           | 136               | 95.1             |
| 4           | 133               | 95.2             |
| 5           | 134               | 95.4             |
| 6           | 133               | 95.5             |
| 7           | 133               | 95.9             |
| 8           | 133               | 95.8             |
| 9           | 130               | 95.7             |
| 10          | 125               | 95.6             |
| 11          | 125               | 95.6             |

|    |     |      |
|----|-----|------|
| 12 | 122 | 95.7 |
| 13 | 119 | 95.8 |
| 14 | 116 | 95.8 |

## Precision

Precision of the System was evaluated by comparing the results from two separate Sensors worn on the same subject at the same time. [Table 27](#) provides data from 146 subjects in the Adult study; [Table 28](#) provides data from 137 subjects in the Paediatric study. For adults, the paired absolute relative difference (PARD) between the two Sensors was 8.1% with coefficient of variation (CV) of 5.8%. For children ages 4 – 5, PARD was 6.7% with CV of 4.8%. For children ages 6 – 17, PARD was 8.3% with CV of 5.9%. Paired absolute difference (PAD) is a measurement of average absolute differences (in mmol/L) between paired GM readings, while PARD is the average of absolute relative differences (in %) between paired GM readings.

**Table 27.** Overall between Sensor Precision (Adult; n=146)

|                 | Mean Coefficient of Variation (%) | Paired Absolute Difference (mmol/L) | Paired Absolute Relative Difference (%) | Number of Paired Readings |
|-----------------|-----------------------------------|-------------------------------------|---|---------------------------|
| Adults ages 18+ | 5.8                               | 0.7                                 | 8.1                                     | 26891                     |

**Table 28.** Overall between Sensor Precision (Paediatric; n=137)

|                      | Mean Coefficient of Variation (%) | Paired Absolute Difference (mmol/L) | Paired Absolute Relative Difference (%) | Number of Paired Readings |
|----------------------|-----------------------------------|-------------------------------------|---|---------------------------|
| Children ages 4 – 5  | 4.8                               | 0.6                                 | 6.7                                     | 248                       |
| Children ages 6 – 17 | 5.9                               | 0.7                                 | 8.3                                     | 10643                     |

## Adverse Events

No device related serious adverse events occurred during the studies. In the Adult study, mild skin irritations (such as erythema, bruising, bleeding and scabbing) and mild pain were reported around the insertion site and adhesive area by a small number of subjects (10 out of 146 or 6.8%).

In the Paediatric study, there were 8 instances of erythema (4 'well-defined redness', and 4 'slight pink'), 5 instances of oedema (3 slight oedema, 2 instances of slight oedema with defined edges), 2 instances of mild bleeding, one instance of mild induration and one instance of mild rash.

## Ascorbic Acid Study

A clinical study was conducted to evaluate the effect of ascorbic acid supplements on Sensor performance. Data from 57 adult subjects with diabetes was collected over a 13-hour period.

Each subject had a one hour baseline phase where venous blood was collected every 10 minutes. After this first hour, a 1000 mg ascorbic acid supplement was given with a meal and venous samples were collected every 20 minutes for the next four hours. A maximum average Sensor bias of 0.5 mmol/L was observed around 3 hours after the 1000 mg ascorbic acid dose. Subjects then received a second 1000 mg ascorbic acid supplement with a meal and the same process was continued for another 4 hours. A third 1000 mg ascorbic acid supplement was then given and study subjects were followed for 4 more hours. After the second dose of ascorbic acid the maximum average Sensor bias increased, with minimal change in Sensor bias after the third dose, suggesting that saturation had occurred by the second 1000 mg dose of ascorbic acid. The maximum average Sensor bias after the three 1000 mg doses of ascorbic acid was less than 1.1 mmol/L.

Bench studies were conducted to test higher levels of ascorbic acid. These studies showed that ascorbic acid doses greater than 1000 mg every four hours could significantly raise Sensor glucose readings.

Customer service: [www.FreeStyleLibre.com](http://www.FreeStyleLibre.com)

Patent: [www.abbott.com/patents](http://www.abbott.com/patents)

FreeStyle, Libre, and related brand marks are marks of Abbott. Other trademarks are the property of their respective owners.

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©2023 Abbott